

## COVID-19 Antigen Test Cassette Product Insert

### INTENDED USE

The SureScreen Diagnostics COVID-19 Antigen Rapid Test is an in vitro immunochromatographic assay for the direct and qualitative detection of SARS-CoV-2 viral nucleoprotein antigens from nasopharyngeal secretions and oropharyngeal secretions from individuals suspected of COVID-19 within the first two weeks of symptom onset. This test is intended for professional use only.

### PRINCIPLE

The SureScreen Diagnostics COVID-19 Antigen Rapid Test detects SARS-CoV-2 viral antigens through visual interpretation of color development. Anti-SARS-CoV-2 antibodies are immobilized on the test region of the nitrocellulose membrane. Anti-SARS-CoV-2 antibodies conjugated to colored particles are immobilized on the conjugated pad. A sample is added to the extraction buffer which is optimized to release the SARS-CoV-2 antigens from specimen. During testing, the extracted antigens bind to anti-SARS-CoV-2 antibodies conjugated to colored particles. As the specimen migrates along the strip by capillary action and interacts with reagents on the membrane, the complex will be captured by the anti-SARS-CoV-2 antibodies at the test region. Excess colored particles are captured at the internal control zone. The presence of a colored band in the test region indicates a positive result for the SARS-CoV-2 viral antigens, while its absence indicates a negative result. A colored band at the control region serves as a procedural control, indicating that the proper volume of specimen has been added and membrane wicking is working.

### MATERIALS

#### Materials Provided

- Individually packed test devices
- Extraction buffer
- Extraction tube
- Nozzle with filter
- Individually packed swabs
- Tube stand
- Package insert

#### Materials Required but Not provided

- Clock, timer or stopwatch

### PRECAUTIONS

- For *in vitro* Diagnostic Use Only.
- Read the Package Insert prior to use. Directions should be read and followed carefully.
- Do not use kit or components beyond the expiration date.
- The device contains material of animal origin and should be handled as a potential biohazard. Do not use if pouch is damaged or open.
- Test devices are packaged in foil pouches that exclude moisture during storage. Inspect each foil pouch before opening. Do not use devices that have holes in the foil or where the pouch has not been completely sealed. Erroneous result may occur if test reagents or components are improperly stored.
- Do not use the Extraction Buffer if it is discolored or turbid. Discoloration or turbidity may be a sign of microbial contamination.
- All patient specimens should be handled and discarded as if they are biologically hazardous. All specimens must be mixed thoroughly before testing to ensure a representative sample prior to testing.
- Failure to bring specimens and reagents to room temperature before testing may decrease assay sensitivity. Inaccurate or inappropriate specimen collection, storage, and transport may yield false negative test results.
- Avoid skin contact with buffer.
- If infection with SARS-CoV-2 is suspected based on current clinical and epidemiological screening criteria recommended by public health authorities, specimens should be collected with appropriate infection control precautions and sent to state or local health departments for testing.
- Viral isolation in cell culture and initial characterization of viral agents recovered in cultures of SARS-CoV-2 specimens are NOT recommended, except in a BSL3 laboratory using BSL3 work practices.

### STORAGE AND STABILITY

- Store the SureScreen Diagnostics COVID-19 Antigen Test Cassette at 2–30°C when not in use.
- **DO NOT FREEZE.**
- Kit contents are stable until the expiration dates marked on their outer packaging and containers.

### SPECIMEN COLLECTION AND STORAGE

#### -Nasopharyngeal swab (NP swab):

- 1) Remove the swab from its packing
- 2) Insert the swab into the nostril parallel to the palate. Rotating against the nasal wall.(to ensure swab contains cells as well as mucus)
- 3) Process the swab as soon as possible after collecting the specimen

#### -Oropharyngeal swab (OP swab):

- 1) Remove the swab from its packing
- 2) Insert the swab completely from the mouth into the throat, centering on the red part of the throat wall and maxillary tonsils and rub the bilateral throat tonsils and throat wall moderately. Avoid touching the tongue and remove the swab
- 3) Process the swab as soon as possible after collecting the specimen

#### Note:


1. Use only synthetic fiber swabs with plastic shafts. Do not use calcium alginate swabs or swabs with wooden shafts, as they may contain substances that inactivate some viruses and inhibit further testing.
2. Swab specimens should be tested as soon as possible after collection. Use freshly collected specimens for best test performance.
3. If not tested immediately, swab specimens may be stored at 2–8°C for 24 hours after collection.
4. Do not use specimens that are obviously contaminate with blood, as it may interfere with the flow of sample with the interpretation of test results.

### TEST PROCEDURE


Bring devices, reagents, and specimens and/or controls to room temperature (15–30°C) before use.

1. For each specimen, open the foil pouch just before testing and remove the test device, and put it on a clean, level surface. Label the tube with the patient identification. For best results, the assay should be performed within one hour.
2. Gently mix extraction buffer. Add 10 drops into the extraction tube.
3. Insert the swab into the extraction tube. Mix well and squeeze the swab 10–15 times by compressing the walls of the tube against the swab. Stand for 2 minutes.
4. Roll the swab head against the inner wall of the tube as you remove it. Try to release as much liquid as possible. Dispose of the used swab in accordance with your biohazard waste disposal protocol.
5. Insert nozzle into sample extraction tube. Invert the tube and add 2 drops of solution into the sample well by gently squeezing the tube.
6. Read results at 15 minutes.

Add 10 drops of buffer

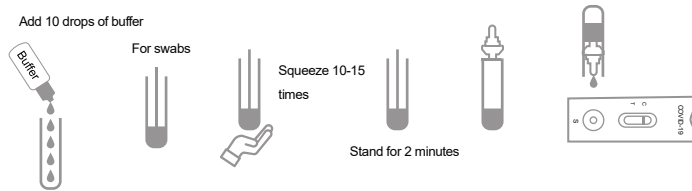

### RESULT INTERPRETATION

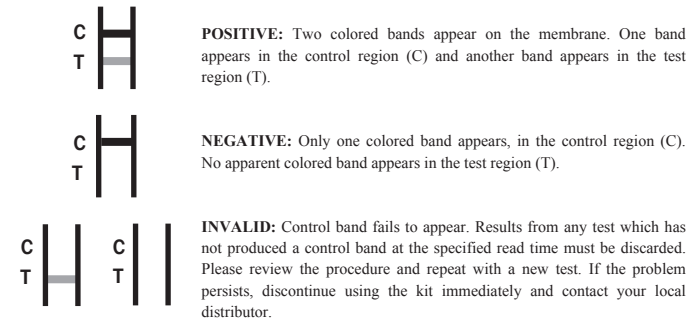

#### Note:

1. The color intensity in the test region (T) may vary depending on the concentration of analytes present in the specimen. Therefore, any shade of color in the test region should be considered positive. Note that this is a qualitative test only and cannot determine the concentration of analytes in the specimen.
2. Insufficient specimen volume, incorrect operating procedure or expired tests are the most likely reasons for control band failure.

### QUALITY CONTROL

#### Internal Procedural Controls

The SureScreen Diagnostics COVID-19 Antigen Rapid Test has built-in (procedural) controls. Each test device has an internal standard zone to ensure proper sample flow. The user should confirm that the colored band located at the “C” region is present before reading the result.

#### External Positive and Negative Controls

Good laboratory practice suggests testing positive and negative external controls to ensure that the test reagents are working and that the test is correctly performed.

### LIMITATIONS OF THE TEST

1. The SureScreen Diagnostics COVID-19 Antigen Rapid Test is for professional in vitro diagnostic use and should only be used for the qualitative detection of SARS-CoV-2 antigen. The intensity of color in a positive band should not be evaluated as “quantitative or semi-quantitative”.
2. Both viable and nonviable SARS-CoV-2 viruses are detectable with the SureScreen Diagnostics COVID-19 Antigen Rapid Test.
3. As with all diagnostic tests, a definitive clinical diagnosis should not be based on the results of a single test but should only be made by the physician after all clinical and laboratory findings have been evaluated.
4. Failure to follow the TEST PROCEDURE and RESULT INTERPRETATION may adversely affect test performance and/or invalidate the test result.
5. Results obtained with this assay, particularly in the case of weak test lines that are difficult to interpret, should be used in conjunction with other clinical information available to the physician.
6. Negative results do not preclude SARS-CoV-2 infection and should be confirmed via molecular assay.

### PERFORMANCE CHARACTERISTICS

#### Analytical Sensitivity (Limit of Detection):

The limit of detection was determined with a quantified SARS-CoV-2 virus and has been evaluated at 2x102.4 TCID50/ml. The limit of detection was also determined with recombinant SARS-CoV-2 nucleoprotein and has been evaluated at 0.4 ng/ml.

#### Clinical Evaluation:

Clinical evaluation was performed to compare the results obtained by the SureScreen Diagnostics COVID-19 Antigen Rapid Test and a real-time Polymerase Chain Reaction (RT-PCR) comparator assay. Patients who presented with symptoms within 14 days were included in the study. The performance was evaluated for 85 prospective clinical specimens which included 43 specimens of nasopharyngeal swab and 42 specimens of oropharyngeal swab. 55 positive specimens and 30 negative specimens were confirmed by RT-PCR. Another 131 retrospective clinical specimens who are either asymptomatic or with mild symptoms (within 14 days of symptom onset) were collected to confirm the specificity. All of these are confirmed negative through PCR. 59 specimens were of nasopharyngeal swab and 72 specimens of oropharyngeal swab. The performance of the SureScreen Diagnostics COVID-19 Antigen Rapid Test Device results based on the various parameters is summarized below:

Table 1: Oropharyngeal Swab Specimen vs. RT-PCR Positive

Days from onset of symptoms	PCR positive at any time	Antigen Positive	PPA	95% Confidence Interval
≤7	6	6	100.00%	51.68% to 100%
8-14	24	22	91.67%	71.52% to 98.54%
≤14	30	28	93.33%	76.49% to 98.83%

Table 2: Nasopharyngeal Swab Specimen vs. RT-PCR Positive

Days from onset of symptoms	PCR positive at any time	SureScreen Diagnostics COVID-19 Antigen Rapid Test		
		Antigen Positive	PPA	95% Confidence Interval
≤7	4	4	100.00%	39.57% to 100%
8-14	21	20	95.24%	74.13% to 99.75%
≤14	25	24	96.00%	77.67% to 99.79%

Table 3: Oropharyngeal and Nasopharyngeal Swab Specimen vs. RT-PCR Positive

Days from onset of symptoms	PCR positive at any time	SureScreen Diagnostics COVID-19 Antigen Rapid Test		
		Antigen Positive	PPA	95% Confidence Interval
≤7	10	10	100%	65.55% to 100%
8-14	45	42	93.33%	80.69% to 98.26%
≤14	55	52	94.55%	83.93% to 98.58%

Table 4: Summary of SureScreen Diagnostic COVID-19 Antigen Rapid Test Clinical Evaluation:

SureScreen Diagnostics COVID-19 Antigen Rapid Test		RT-PCR		Total
		Positive	Negative	
	Positive	52	0	52
	Negative	3	161	164
	Total	55	161	216

Relative Sensitivity: 94.55% (83.93%–98.58%)\*

Relative Specificity: 100% (97.1%–100%)\*

Overall Agreement: 98.61% (95.66%–99.64%)\*

\*95% Confidence Interval

#### Cross Reactivity:

Cross reactivity with the following organisms has been studied. Samples positive for the following organisms were found negative when tested with the SureScreen Diagnostic COVID-19 Antigen Rapid Test (Nasopharyngeal/Oropharyngeal Swab).









HCoV-HKU1	Influenza A (H5N1)	Coxsackie virus A16
HCoV-OC43	Influenza A (H7N9)	Norovirus
HCoV-NL63	Influenza A (H7N7)	Mump virus
HCoV-229E	Influenza B Victoria lineage	<i>Legionella pneumophila</i>
Measles virus	Influenza B Yamagata lineage	<i>Mycoplasma pneumoniae</i>
Streptococcus pneumoniae	Respiratory syncytial virus	<i>Chlamydia pneumoniae</i>
Epstein-Barr virus	Adenovirus	<i>Streptococcus pyogenes</i>
Bordetella parapertussis	Parainfluenza 1/2/3 virus	<i>Streptococcus agalactiae</i>
InfluenzaA (H1N1) pdm09	Human metapneumovirus	Group C <i>Streptococcus</i>
Influenza A (H3N2)	Rhinovirus	<i>Staphylococcus aureus</i>

#### Interfering Substances

The following substances, naturally present in respiratory specimens or that may be artificially introduced into the respiratory tract, were evaluated at the concentrations listed below. None of them were found to affect test performance of the SureScreen Diagnostics COVID-19 Antigen Rapid Test Device.

Substance	Concentration	Substance	Concentration
3 OTC nasal sprays	10%	Guaiacol glyceryl ether	20 mg/ml
3 OTC mouthwashes	10%	Mucin	1%
3 OTC throat drops	10%	Mupirocin	250 µg/ml
4-acetamidophenol	10 mg/ml	Oxymetazoline	10 mg/ml
Acetylsalicylic acid	20 mg/ml	Phenylephrine	10 mg/ml
Albuterol	20 mg/ml	Phenylpropanolamine	20 mg/ml
Chlorpheniramine	5 mg/ml	Relenza <sup>®</sup> (zanamivir)	20 mg/ml
Dexamethasone	5 mg/ml	Rimantadine	500 ng/ml
Dextromethorphan	10 mg/ml	Tamiflu <sup>®</sup> (oseltamivir)	100 mg/ml
Diphenhydramine	5 mg/ml	Tobramycin	40 mg/ml
Doxylaminesuccinate	1 mg/ml	Triamcinolone	14 mg/ml
Flunisolide	3 mg/ml		

#### GLOSSARY OF SYMBOLS

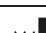 Consult instructions for use	 Test per Kit	 Catalogue #
 Store between 2°C to 30°C	 Use by	 Do Not Reuse
 Manufacturer	 Lot Number	



Number: 1110032811  
Effective Date: 04/09/2020

#### LITERATURE REFERENCES

- Forni, D., Cagliani, R., Clerici, M. & Sironi, M. Molecular evolution of human coronavirus genomes. Trends Microbiol. 25, 35–48 (2017).
- Ithete, N. L. et al. Close relative of human Middle East respiratory syndrome coronavirus in bat, South Africa. Emerg. Infect. Dis. 19, 1697–1699 (2013).



**Wolflabs**

# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

