



Product Instructions

RS96010BPW / RS96010DE / RS96010LET /

RS9605LET / RS9604LET / RS9601LET / RS96010NB / RS9604NB

-  **(EN) Swab-Sampler**
-  **(FR) Ecouvillon**
-  **(DE) Probestupfer**
-  **(IT) Tampone per campionamento**
-  **(ES) Hisopo muestrador**
-  **(NL) Swab-Sampler**
-  **(SV) Swab-Sampler**
-  **(DA) Svaber til prøveudtagning**
-  **(NO) Prøvesvaber**
-  **(FI) Swab-Sampler Näytteenkeräystikku**
-  **(PT) Swab-Sampler**
-  **(EL) Στυλεός Δειγματοληψίας**
-  **(PL) Próbnik do wymazów**
-  **(HU) Swab-Sampler kenetvevő**
-  **(CZ) Stěrový vzorkovač**
-  **(RO) Tampon de recoltare**
-  **(RU) Ультратампон**
-  **(TR) Svab Örnekleyici**
-  **(JA) スwabサンプルラー**
-  **(ZH) 采样棒**
-  **(TH) ก้านสวอป**
-  **(KO) 스왑-샘플러**

Product Instructions

Swab-Sampler

PRODUCT DESCRIPTION

Pre-measured volume of enrichment broth or neutralizing solution in a tube with a screw-cap and attached swab.

For Laboratory Use Only. Material Safety Data Sheet available.

SAFETY

The user should read, understand, and follow all safety information in the instructions. Retain the safety instructions for future reference.

⚠ WARNING Indicates a hazardous situation, which, if not avoided, could result in death or serious injury and/or property damage.

NOTICE Indicates a potentially hazardous situation which, if not avoided, could result in property damage.

⚠ WARNING

To reduce the risks associated with environmental contamination:

- Follow current industry standards and local regulations for disposal of contaminated waste.

To reduce the risks associated with exposure to biohazards:

- Dispose of samples according to all applicable government regulations and applicable laboratory procedures.

To reduce the risk associated with false negatives resulting in the use of contaminated carcass or environmental surfaces for food or beverage products:

- Always reference package label for storage instruction and expiration date.
- Always reference product instruction for usage.

To reduce the risks associated with exposure to chemicals and biohazards:

- Always follow standard good laboratory safety practices (GLP¹ or ISO 17025²), including proper containment procedures, wearing appropriate protective apparel while handling testing materials and test samples.

NOTICE

To reduce the risk of false-positive results due to cross contaminated carcass or environmental surfaces for food or beverage products that may result in re-testing or the rejection of food or beverage product:

- Do not touch the swab or swab stick to any unintended surface.

To reduce the risk of cross contamination from reuse of sample handling device:

- Do not use the same swab more than once.

USER RESPONSIBILITY

Users are responsible for familiarizing themselves with product instructions and information. Visit our website at www.3M.com/foodsafety, or contact your local 3M representative or distributor for more information.

When selecting a test method, it is important to recognize that external factors such as sampling methods, testing protocols, sample preparation, handling, and laboratory technique may influence results.

It is the user's responsibility in selecting any test method or product to evaluate a sufficient number of samples with the appropriate matrices and microbial challenges to satisfy the user that the chosen test method meets the user's criteria.

It is also the user's responsibility to determine that any test methods and results meet its customers' and suppliers' requirements.

As with any test method, results obtained from use of any 3M Food Safety product do not constitute a guarantee of the quality of the matrices or processes tested.

LIMITATION OF WARRANTIES / LIMITED REMEDY

EXCEPT AS EXPRESSLY STATED IN A LIMITED WARRANTY SECTION OF INDIVIDUAL PRODUCT PACKAGING, 3M DISCLAIMS ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. If any 3M Food Safety Product is defective, 3M or its authorized distributor will, at its option, replace or refund the purchase price of the product. These are your exclusive remedies. You must promptly notify 3M within sixty days of discovery of any suspected defects in a product and return it to 3M. Please call Customer Service (1-800-328-1671 in the U.S.) or your official 3M Food Safety representative for a Returned Goods Authorization.



LIMITATION OF 3M LIABILITY

3M WILL NOT BE LIABLE FOR ANY LOSS OR DAMAGES, WHETHER DIRECT, INDIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES, INCLUDING BUT NOT LIMITED TO LOST PROFITS. In no event shall 3M's liability under any legal theory exceed the purchase price of the product alleged to be defective.

STORAGE AND DISPOSAL

Storage: Refer to package label for storage information.

Disposal: Follow current industry standards and local regulations for disposal of contaminated waste. Consult the Material Safety Data Sheet for additional information.

Product Description	Volume	REF
3M™ Swab-Sampler with Buffered Peptone Water Broth	10 mL	RS96010BPW
3M™ Swab-Sampler with D/E Neutralizing Broth	10 mL	RS96010DE
3M™ Swab-Sampler with Lethen Broth	10 mL	RS96010LET
3M™ Swab-Sampler with Lethen Broth	5 mL	RS9605LET
3M™ Swab-Sampler with Lethen Broth	4 mL	RS9604LET
3M™ Swab-Sampler with Lethen Broth	1 mL	RS9601LET
3M™ Swab-Sampler with Neutralizing Buffer	10 mL	RS96010NB
3M™ Swab-Sampler with Neutralizing Buffer	4 mL	RS9604NB

INSTRUCTIONS FOR USE

1. Unscrew cap from tube.
2. Aseptically remove swab from tube.
3. Aseptically swab across the entire sampling surface^{3,4,5,6} while rotating swab.
4. Return swab to tube.
5. Repeat step 2. Change direction 90° and aseptically swab the same sampled surface while rotating swab. Repeat step 4.
6. Repeat step 2. Change direction 45° and aseptically swab the same sampled surface while rotating swab. Repeat step 4.
7. Screw cap tight to close.
8. Following user established procedures, remove any remaining enrichment broth or neutralizing solution residue from the sampled surface.

REFERENCES

1. U.S. Food and Drug Administration. Code of Federal Regulations, Title 21, Part 58. Good Laboratory Practice for Nonclinical Laboratory Studies.
2. ISO/IEC 17025. General requirements for the competence of testing and calibration laboratories.
3. American Public Health Association - Compendium of Methods for the Microbiological Examination of Foods - Chapter 3: Microbiological Monitoring of the Food Processing Environment, 4th edition.
4. US Food and Drug Administration - Bacteriological Analytical Method (available online at <http://www.fda.gov/Food/ScienceResearch/LaboratoryMethods/BacteriologicalAnalyticalManualBAM/default.html>).
5. United States Department of Agriculture - Microbiological Lab Guidebook - Chapters 4.04, 5.04, and 8.07.
6. American Public Health Association - Standard Methods for the Examination of Dairy Products - Chapter 13: Microbiological Tests for Equipment, Containers, Water, and Air, 17th edition.

Refer to the current version of the standard methods listed above.

EXPLANATION OF SYMBOLS












Consult instructions for use



The lot in a box and the hourglass symbols are symbols that represent lot number and expiration date. The hourglass is followed by a year, month, and day which represents the expiration date (year, month, and day: YYYY- MM-DD). The entire line after the LOT represents the lot number (YYYY- MM AZ).



Storage conditions

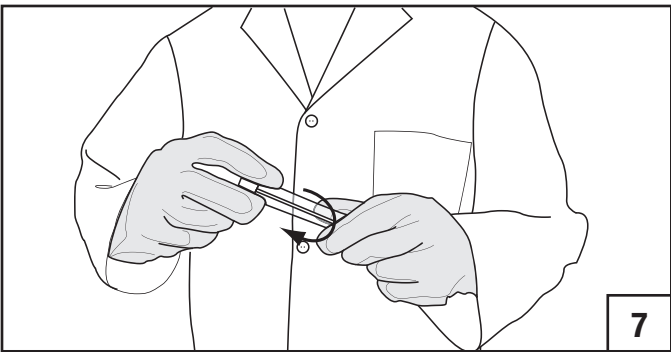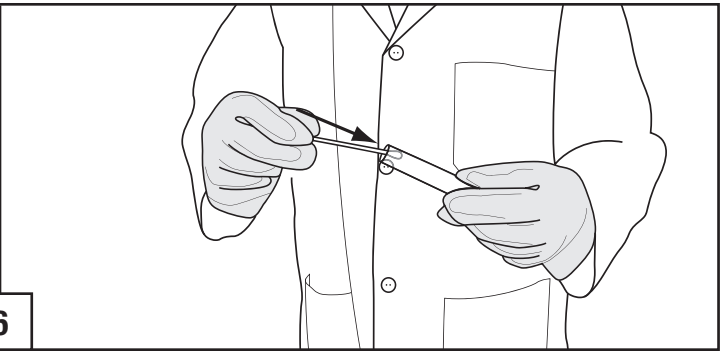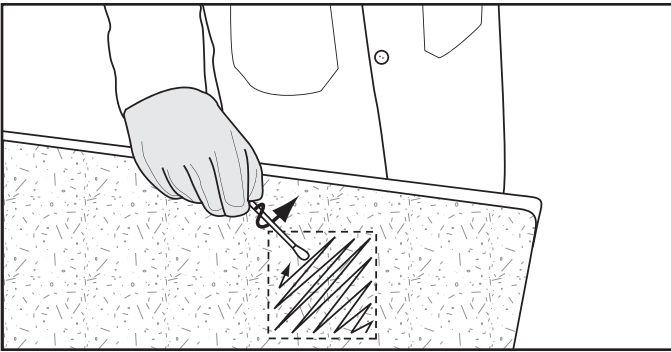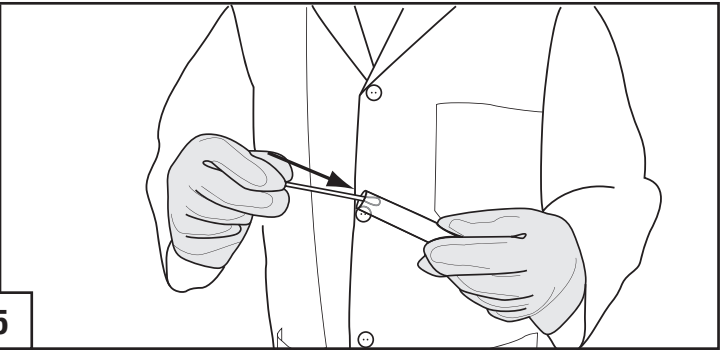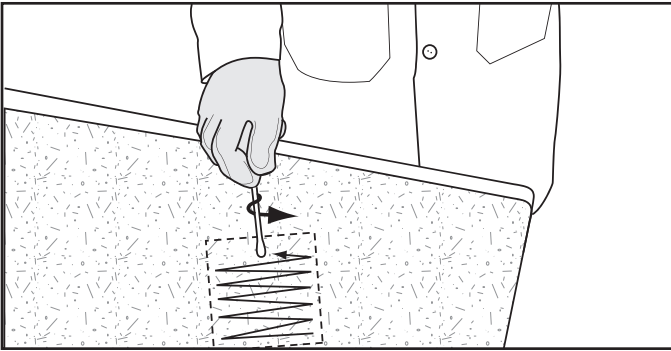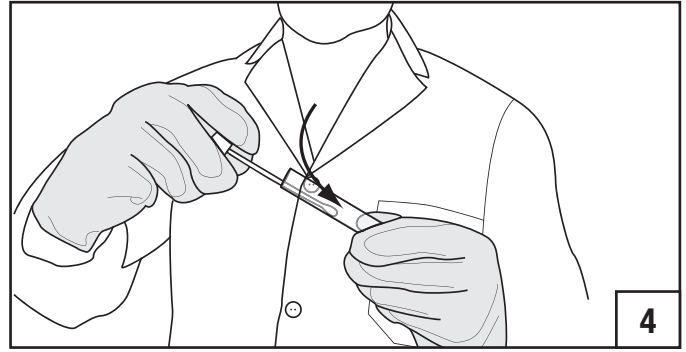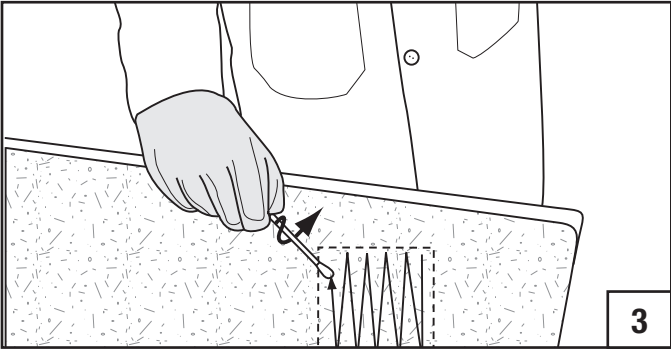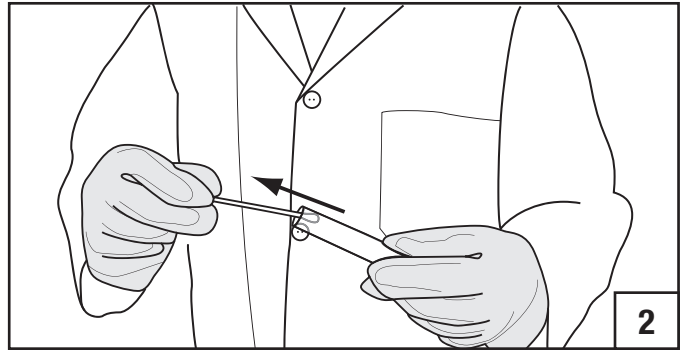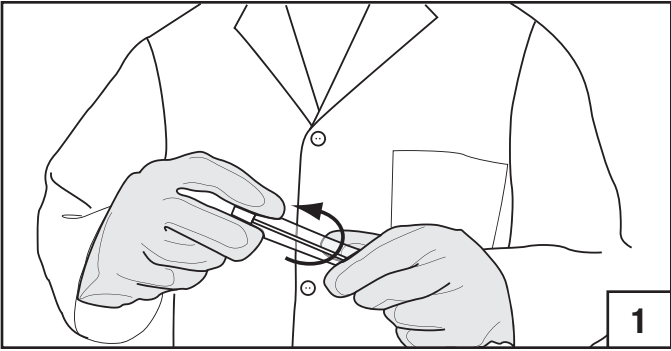



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.