




Built to withstand the harshest working conditions, Aqua washdown scales are rated IP67 and feature a large removable stainless steel weighing pan. Ideal for food production, catering and industrial environments, precision load technology enables improved performance and greater efficiency during weighing and counting tasks. The dual front and rear displays enhance visibility and productivity for the user. Checkweighing is simple, with separate high, ok and low coloured LED indicators that change colour to alert the operator if a sample is under, over or within preset limits. Power is supplied by AC adapter or rechargeable battery, and auto power-off conserves battery life.

*Get the Aqua wet without worry, the IP67 scale performs effortlessly and accurately no matter how rigorous the application*

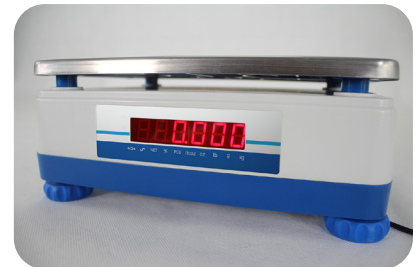

## Features

- IP67-rated for protection against water and dust
- External calibration allows for verification and adjustment with weights
- Sealed keypad protects against dirt and spills
- Overload protection helps prevent damage to internal components
- Solid ABS housing is durable for industrial use
- Heavy-duty rubber feet ensure stability of unit
- Two LED displays for dual use, with 0.8" / 20mm-high (front) and 0.6" / 15mm-high (rear) digits
- Checkweighing LEDs clearly show under, over and within acceptable limits
- Parts counting optimisation automatically refines piece weight as parts are added
- Zero-tracking feature ensures display returns to zero reading
- The auto power-off feature conserves battery life
- Rechargeable battery included for operation almost anywhere
- AC adapter included

## Applications

- Weighing
- Parts counting
- Percentage weighing
- Checkweighing

### About ADAM:

For more than 50 years, professionals worldwide have relied on Adam Equipment for an extensive selection of dependable and affordable weighing equipment. Adam provides the right balance of speed, performance and value.

## Product Details

Model	ABW 4	ABW 8	ABW 16	ABW 32
Capacity	4kg	8kg	16kg	32kg
Readability	0.5g	1g	2g	5g
Repeatability (S.D.)	1g	2g	5g	10g
Linearity (+/-)	1g	2g	5g	10g
Pan Size	245 x 180mm			300 x 210mm
Weighing Units	g, kg, lb, oz, lb:oz			
Stabilisation Time (sec)	1			
Calibration	External			
Display	LED with 20mm-high (front) and 15mm-high (rear) digits			
Power Supply	12VDC 800mA adapter, Internal rechargeable battery			
Operating Temp.	0° to 40°C			
Housing	ABS Plastic			
Overall Dim.	255 x 290 x 105mm			300 x 310 x 105mm

## Accessories

### Item Number.....Description

- 3022014061..... In-use cover
- 2020014062 ..... In-use cover (pack of 5)
- 2020014063 ..... In-use cover (pack of 10)
- 3002014371 ..... Hard carrying case
- 303149759..... Large Scoop (complete with fitting to scales)
- 303149760..... Small Scoop (complete with fitting to scales)
- 700660290..... Calibration certificate

## Key Accessories



### Scoops

Metal scoops are designed to make portioning and food measurement easy. Two sizes available: small and large.



### Hard Carrying Case

Allows the Aqua to be placed inside to protect it during transportation. The hard shell case is durable and comes with a lock so that you can secure the scale in the case as needed. Die cut foam insert provides room for the scale, power supply and top pan.

Supplied by



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.