



SL-CLASSIC-UB ◀ 059 cm ▶

NEW

5 Years
Warranty



Safe and approved storage of aggressive, non flammable hazardous materials in working areas

Function / construction:

- **robust construction and longevity:** made out of highly resistant melamin resin special plates
- **user-friendly:** access to all containers, 135° opening of the wing door, interior of the cabinet completely visible
- **no unauthorised use:** door lockable with cylinder lock and locking state indicator (red/green)
- **mobile:** castors with plinth (optionally available), easy moving of the cabinet
- **ventilation:** integrated air ducts ready for connection (DN 50) to a technical exhaust system

Available equipment

- pull-out shelves with removable plastic sump

Cabinet for acids and alkalis

SL-CLASSIC-UB model SL.080.110.UB.2T

in Light grey (RAL 7035), interior equipment with 4 x pull-out shelf (panel/PP)

Order No. 33200-001-30621



castors with plinth

stable – easy moving of the cabinet

safe fixing of the cabinet – integrated locking mechanism

flexible – height adjustment to various laboratory furniture



Cabinet for acids and alkalis

SL-CLASSIC-UB model SL.080.059.UB.TR

in Light grey (RAL 7035), interior equipment with 2 x pull-out shelf (panel/PP)

Order No. 33202-001-30625

Technical data

External dimensions W x D x H

mm 590 x 570 x 800

1100 x 570 x 800

Internal dimensions W x D x H

mm 490 x 530 x 760

all compartments: 1000 x 530 x 760

Weight without interior equip.

kg 48

73

Extraction air

DN 50

50

Air change **30 times**

m³/h 6

9

Diff. pressure (air change **30 times**)

Pa 5

10

Distributed load

kg/m² 157

98


Maximum load

kg 100


100




SL.080.110.UB.2T Front view

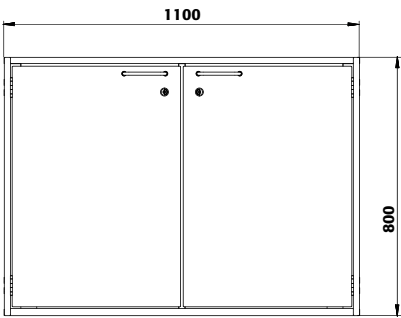


Side view

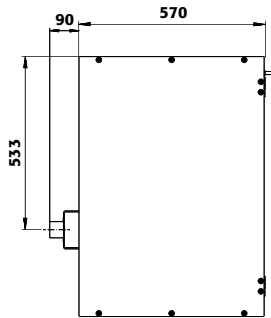


Top view

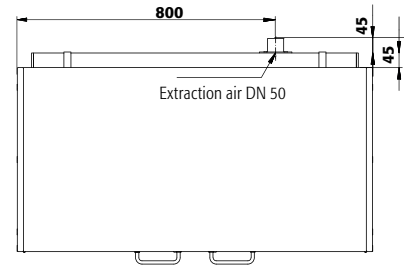


SL.080.059.UB.TR Front view

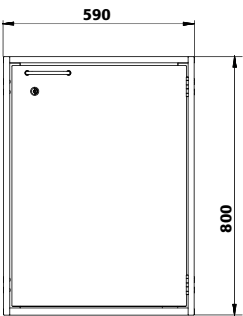

Side view

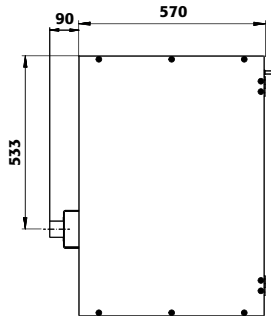

Top view

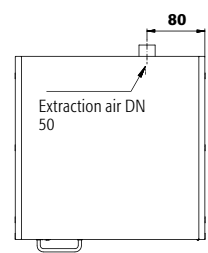

SL-LINE

Order example - how to order correctly



Model	Information / equipment	Order No.
SL.080.110.UB.2T	without interior equipment, Underbench cabinet, with 2 wing doors	33200
SL.080.059.UB.TR	without interior equipment, Underbench cabinet, right hinged, with wing door	33202

Body colour	Door colour	Order No.
Light grey RAL 7035	Light grey RAL 7035	001 4

Interior equipment packages (can only be ordered in combination with the cabinet)	Material	Order No.
Width 1100 mm		
left: 1x pull-out shelf (V=11.0L)	panel/PP	30620 2
right: 1x pull-out shelf (V=11.0L)		
left: 2x pull-out shelf (V=11.0L)	panel/PP	30621 2
right: 2x pull-out shelf (V=11.0L)		
Width 590 mm		
1x pull-out shelf (V=11.0L)	panel/PP	30624 2
2x pull-out shelf (V=11.0L)	panel/PP	30625 2

Interior equipment	Material	W x D x H (mm)	1)	2)	3)	Order No.
pull-out shelf, with PP sump	panel/PP grey	590, 1100 and 1400	25.00	18.00	16.36	30770 2
castors with plinth, height = 30 mm	sheet steel powder-coated RAL 7035	1100				30362 2
castors with plinth, height = 30 mm	sheet steel powder-coated RAL 7035	590				30478 2

Accessories	Order No.
extraction module without exhaust air monitoring for wall mounting	17178 2

* usable storage space W x D (mm)
 1) load capacity (kg), 2) retention capacity (ltr), 3) maximum container volume in the cabinet (ltr)



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.