

Benchmark 

myBLOCK™ Digital Dry Baths

- + **Digital temperature control set and walk away**
- + **Hinged lid improves temperature uniformity and conserves energy**
- + **External temperature probe (optional) displays and controls actual sample temperature**
- + **Available with Quick-Flip™ universal blocks for microtubes (0.2 to 2.0ml), PCR strips & PCR plates**

Flip Over



myBLOCK I
Single Block Capacity



myBLOCK II
Dual Block Capacity





myBlock I (single chamber dry bath)

Benchmark's myBlock series is a "next generation" product line that introduces several major advantages over traditional dry baths. They are the first digital dry baths to include advanced microprocessor controls, timed or continuous operation and a removable hinged lid. In addition, an optional external temperature feedback probe can be positioned directly inside a sample tube, providing real-time control and monitoring of **actual sample temperature**.

The myBlock's transparent lid helps to promote better uniformity and minimize temperature stratification in the block by maintaining a warm air jacket layer. It also helps to reduce energy consumption by minimizing heat loss to the environment. The hinged lid can be flipped open (to the rear), or easily removed (lifted off) to accommodate taller tubes.

Both the single and dual block models can be ordered in a configuration that includes our unique Quick-Flip™ blocks. These universal blocks accept tubes from 0.5 to 2.0ml and, with just a quick flip to the other side, accommodate 0.2ml tubes, PCR strips or (in the dual block model only) a PCR plate. Other blocks are available for tubes up to 50ml.

Technical Data:

Temperature Range:	Ambient +5 to 105°C
Temp. Accuracy / Uniformity:	± 0.3°C / ± 0.1°C
Temperature Increments:	0.1°C
Timer:	1 min to 99 hr 59 min.
Dimensions (w x d x h):	
Single:	7 x 10 x 5 in. / 18 x 25 x 13 cm
Dual:	6.3 x 14 x 5 in. / 16 x 36.5 x 13 cm
Weight:	10 lbs. / 4.5 kg
Electrical:	115V or 230V, 50-60Hz
Warranty:	2 years

available with:

Quick-Flip™
universal micro-block



Universal blocks for:

- ✓ 0.5ml tubes
- ✓ 1.5ml tubes
- ✓ 2.0ml tubes
- ✓ 0.2ml tubes
- ✓ 0.2ml strips
- ✓ PCR Plate

Item: BSH5002-2B


Three possible block configurations:
(Shown in dual block model)

Side 1 & Side 1



48 x 1.5/2.0ml
28 x 0.5ml

Side 1 & Side 2

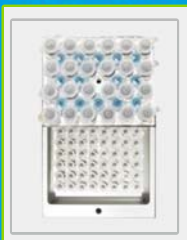

24 x 1.5/2.0ml
14 x 0.5ml
48 x 0.2ml

Side 2 & Side 2



96 x 0.2ml or
One PCR Plate*

*For use in myBlock II only - requires dual block chamber.

Ordering Information

- BSH5001-1B* **myBlock I** Single chamber dry bath with 1 Quick-Flip block (BSWCMB) for tubes (0.2 to 2.0ml) and PCR strips
- BSH5002-2B* **myBlock II** Dual chamber dry bath with 2 Quick-Flip blocks (BSWCMB) for tubes (0.2 to 2.0ml), PCR strips / PCR Plate
- BSH5001* **myBlock I** single chamber dry bath, without block
- BSH5002* **myBlock II** dual chamber dry bath, without blocks
- BSH-TP1 Optional external temperature probe. For monitoring & controlling actual sample temperature. (Probe can be placed in the sample tube with direct feedback to the control circuit and temperature display.)

* Item is 115V. To order in 230V, add -E to catalog number.

Available Block Accesories:

- BSWCMB **Quick-Flip** combination block, Side one holds: 24 x 1.5/2.0ml and 14 x 0.5ml, Side two holds: 6 x PCR strips, 48 x 0.2ml or one PCR plate*
- BSW1500 Block, 24 x 1.5ml
- BSW1520 Block, 24 x 1.5/2.0ml
- BSW05 Block, 24 x 0.5ml
- BSW10 Block, 20 x 10mm†
- BSW13 Block, 20 x 13mm (or 20 x 5/7ml)†
- BSW1516 Block, 12 x 15-16mm (or 12 x 10ml)†
- BSW02 Block, 48 x 0.2ml tubes or 6 PCR strips
- BSW15 Block, 12 x 15ml†
- BSW50 Block, 5 x 50ml†
- BSWPCR1 Block, 1 x 96-well PCR Plate
- BSWPCR2 Block, 1 x 96-well PCR Plate (myBlock II only)
- BSWMT Block, 1 x Micro Titer Plate (myBlock II only)
- BSW01 Block, solid - for slides or custom machining
- BSWHEM 24 x hematocrit tubes†

* Two blocks required for use with PCR plates. † Tall block (3.25"), height prevents use with lid.



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.