

ORBi-SHAKER™ XL

MAXIMUM CAPACITY, BENCHTOP MODEL



CE

with **MAGiC CLAMP™**
PATENTED

- + High Capacity, benchtop model
- + Circular 19mm orbit for aeration and mixing
- + Exchange flask clamps in one second with MAGic Clamp™ platform.

Technical Data

Speed Range:	30 to 300 rpm
Speed Increment:	1 rpm
Orbit:	19mm (3/4")
Max. Capacity:	5 x 2L or 9 x 1L Flasks
Max. Load:	65 lbs (30kg)
Platform Dimensions:	18 x 18 in. / 45 x 45 cm (Included)
Timer:	1 min. to 99 hours / continuous
Weight:	65 lbs / 30 kg
Dimensions:	18.5 x 21.25 x 5 in. / 47 x 54 x 12.5 cm
Electrical:	120V or 230V, 50-60Hz

The Orbi-Shaker XL offers flask capacity that far exceeds most benchtop models. With a large 18 x 18 inch platform and a powerful, beltless motor drive, the Orbi-Shaker XL is capable of mixing up to five 2L flasks or nine 1L flasks simultaneously. The included, non-slip rubber mat platform is ideal for low speed applications, while the convenient MAGic Clamp™ platform is also available for high speed mixing of Erlenmeyer flasks, tube racks or microplates. The MAGic Clamp™ allows for instant flask clamp exchange without the use of tools, screws or any other hardware.

Ordering Information




BT1011* Orbi-Shaker XL, High capacity, benchtop shaker
Includes 18x18" platform with non-slip mat, 115V

Optional Platforms and Accessories:

H1010-MR	MAGic Clamp universal platform, (18x18")
H1000-MR-50	Flask Clamp, 50ml (max. 49)
H1000-MR-125	Flask Clamp, 125ml (max. 36)
H1000-MR-250	Flask Clamp, 250ml (max. 25)
H1000-MR-500	Flask Clamp, 500ml (max. 16)
H1000-MR-1000	Flask Clamp, 1000ml (max. 9)
H1000-MR-2000	Flask Clamp, 2000ml (max. 5)
H1000-MR-1550	Tube Rack, 20x50 & 30x15 ml tubes (max. 4)
H1000-MR-T15	Tube Rack, 32x15 tubes (max. 8)
H1000-MR-T50	Tube Rack, 15x50 tubes (max. 8)
H1000-MR-T600	Bottle Rack, 4x500/600ml conical bottles (max. 4)
H1000-MR-T15H	Tube Rack, 12 x 1.5/2.0ml horizontal (max. 6)
H1000-MR-T150H	Tube Rack, 4 x 15ml horizontal (max.8)
H1000-MR-T500H	Tube Rack, 2 x 50ml horizontal (max.8)

H1010-P-SP Spring loaded platform, for uncommon vessels
H1010-P-MP 6 x microplate platform (standard or deepwell)

*To order item in 230V, add (-E) to item number.

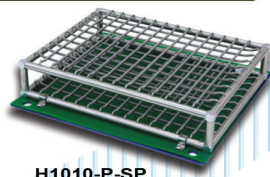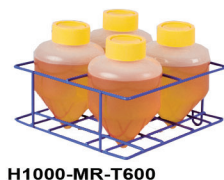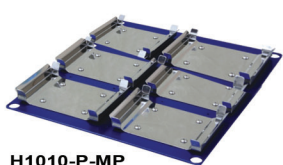



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.