

# Model BD-S 115 | Standard-Incubators with natural convection

The new products in the Solid.Line impress with their tried-and-tested BINDER quality and reliability. They are suited to incubation. The attributes of the Solid.Line are geared towards the application in question and therefore used in research and quality assurance.

## BENEFITS

- Easy handling
- Accurate temperature control
- Simple and ergonomic door opening



Model 115

## MAIN FEATURES

- Temperature range: ambient temperature +5 °C to 70 °C
- APT.line™ preheating chamber technology
- Natural convection
- Adjustable exhaust air flap
- Controller with timer function
- Inner door made of tempered safety glass
- 1 chrome-plated rack, incl. shelf supports
- Class 3.1 integrated independent temperature safety device (DIN 12880) with visual alarm

## ORDERING INFORMATION


Interior volume [L]	Voltage	Option model	Version	Art.-No.
Model BD-S 115				
118	230 V 1~ ph 50/60 Hz	Standard	BDS115-230V	9090-0022
	120 V 1~ ph 60 Hz	Standard	BDS115UL-120V	9090-0023

## TECHNICAL DATA

Description	BDS115-230V	BDS115UL-120V
Article Number	9090-0022	9090-0023
Performance Data Temperature		
Temperature range 5 °C above ambient temperature to [°C]	70	70
Temperature variation at 37 °C [± K]	0.5	0.5
Temperature fluctuation at 37 °C [± K]	0.2	0.2
Recovery time after 30 seconds door open at 37 °C [min]	15	15

Description	BDS115-230V	BDS115UL-120V
Article Number	9090-0022	9090-0023
<b>Electrical data</b>		
Rated Voltage [V]	230	120
Power frequency [Hz]	50/60	60
Nominal power [kW]	0.35	0.35
Unit fuse [A]	6.3	6.3
Phase (Nominal voltage)	1~	1~
<b>Internal Dimensions</b>		
Depth [mm]	390	390
Height [mm]	550	550
Width [mm]	550	550
<b>Housing dimensions not incl. fittings and connections</b>		
Width net [mm]	764	764
Height net [mm]	735	735
Depth net [mm]	616	616
<b>Measures</b>		
Wall clearance back [mm]	160	160
Wall clearance sidewise [mm]	100	100
Interior volume [L]	118	118
Net weight of the unit (empty) [kg]	52.5	52.5
permitted load [kg]	75	75
Load per rack [kg]	20	20
<b>Doors</b>		
Inner doors	1	1
Unit doors	1	1
<b>Environment-specific data</b>		
Energy consumption at 37 °C [Wh/h]	25	25
<b>Fixtures</b>		
Number of shelves (std./max.)	1/5	1/5

### DIMENSIONS incl. fittings and connections [mm]

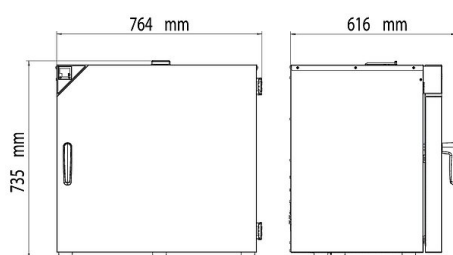

### OPTIONS AND ACCESSORIES

Designation	Description	*	Art.-No.
Data Logger Kit	T 220: For continuous temperature logging from -90 °C to 220 °C. The kit includes 1 data logger, Pt 100 sensor with 2 m extension cable and 1 magnetic fixture for mounting to the BINDER unit	19	8012-0715
Data Logger Software	LOG ANALYZE software kit, configuration and evaluation software for all BINDER Data Logger Kits (incl. USB data cable)	19	8012-0821
Grid shelf	Incl. shelf supports chrome plated	-	8012-1808
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	1002-0016
Rubber pads	set anti-slip feet	-	8012-1887

\* Notes > See last page

## SERVICES

Designation	Description	*	Art.-No.
<b>Installation services</b>			
Installation	and set up of unit at operating location, connect to existing connections	13, 18	DL10-0100
Instruction	unit function instructions for operation and programming of the controller	18	DL10-0500
<b>Preventive maintenance</b>			
Preventive maintenance	Executive of equipment inspection according to maintenance plan	14, 18	DL20-0200
<b>Calibration services</b>			
Calibration temperature incl. certificate	including certificate, one measuring point in center of chamber at specified temperature	14, 16, 17, 18	DL30-0101
Temperature measurement 27 measuring points	including certificate, 27 measuring points at specified temperature	14, 16, 17, 18	DL30-0127
Temperature measurement 9 measuring points	including certificate, 9 measuring points at specified temperature	14, 16, 17, 18	DL30-0109
<b>Warranty service</b>			
Extension of the warranty from 2 to 3 years	beginning with the date of delivery, wearing parts are not included	-	DL01-6041

\* Notes › See last page



**Wolflabs**

# Wolf Laboratories Limited

[www.wolflabs.co.uk](http://www.wolflabs.co.uk)

Tel: 01759 301142

Fax: 01759 301143

[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)



**Use the above details to contact us if this literature doesn't answer all your questions.**

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

