


CAPP Rondo

TUBE ROTATOR

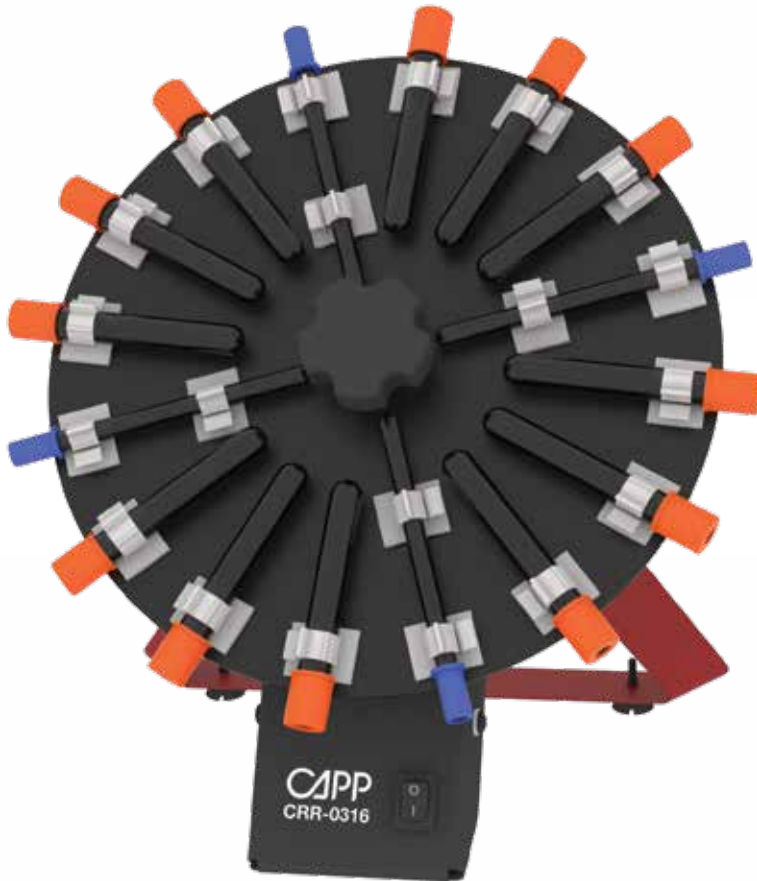

CAPP Rondo Tube Rotator has a strong, metal design which, along with the rubber feet, ensures its total stability during operation.

The compact design uses a minimal amount of benchtop space.

The gentle but effective mixing occurs at a fixed speed of 30 RPM with an accuracy of +/- 2 RPM.

CAPP Rondo Tube Rotator has a maximum capacity of 16 tubes, consisting of 4 ESR and 12 blood collection tubes. Its maximum operational weight capacity is 1.5kg. The disc plate is fixed at an angle of 38 degrees and is easily removed for cleaning or replacement.

CAPP Rondo Tube Rotator operates silently in the continuous mode.

Specifications:

Max Capacity	16 Tubes (4 ESR Tubes and 12 Blood Collection Tubes)
Motor Type	AC Synchronous
Speed	Fixed Speed of 30 RPM
Operation Mode	Continuous / Infinite
Speed Accuracy	+/- 2 RPM
Angle Of Disc Plate	Fixed angle of 38 degrees
Noise Level	Noiseless operation
Load Capacity	1.5kg
Dimensions	278x183x252mm
Protection Class	IP 21

Ordering Information

Cat No.	Description
CRR-0316	CAPPRondo Tube Rotator for 16 blood tubes
CRR-0316-110	CAPPRondo Tube Rotator for 16 blood tubes, 110V





WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.