

Clifton Boiling Baths

NE1B Series Boiling Baths

Clifton Range[®]



NE1B-14




Designed for applications that require continuous boiling. This range of baths can be regulated from gentle simmer to vigorous boil. With the heater directly immersed in water it provides even heat distribution across the tank base. The perforated tray ensures an even spread of heat throughout the chamber. Boiling bath comes fitted with Clear lid.

Comprising of a high quality corrosion resistant stainless steel tank with bright clean finish housed in a durable painted outer case. Painted surfaces feature anti-bacterial finish with good chemical resistance.

The analogue dial has a proportional 1-10 scale, and the energy regulator provides control from gentle simmer to vigorous boil.

	NE1B - 14	NE1B - 18W	NE1B - 28
GENERAL			
Temp. Range	Simmer to boiling at 100°C		
Temp. Control	Analogue dial, energy regulator		
Heater Power (220-240V)	2500W		
Water Level Alarm	Features low liquid level float device, heating automatically suspended		
Voltage	230V / 50-60Hz		
Drain (12mm Bore)	Drain tap fitted		
DIMENSIONS			
Tank Capacity	14 Litre	18 Litre	28 Litre
External (mm)	362w x 332d x 290h		538w x 332d x 290h
Internal (mm)	325w x 300d x 150h	325w x 300d x 200h	500w x 300d x 200h

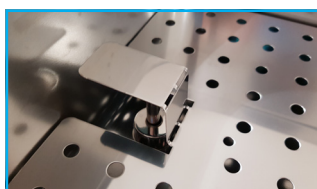

Low liquid float device

A low liquid float device protects the heaters and a red warning light on the front panel alerts users to refill the bath.



NE1B-18W

Raise and lower sterilising tray along with pair of guide rail assemblies comes as standard.

NE1B Series Boiling Baths

Accessories



Stainless Steel Gabled Lid with Insulated Handle

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	N/A	N/A	SL1-22SS



Stainless Steel Flat One Piece Lid with Insulated Handle

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	LD-14	N/A	LD-22



Stainless Steel Hinged Gabled Lid with Insulated Handle

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	SL1-14H/SS	N/A	N/A



Clear Hinged Gable Lid (Fitted as standard)

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	GL1-14		GL1-22



Clear Lift Off Gable Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	GL1-14AF	N/A	GL1-22AF



Stainless Steel Lid with Concentric Rings

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	SLR1-14	N/A	SLR1-22
	SLR2-14		

SLR1-4: 2 x 114mm holes

SLR1-14: 4 x 114mm holes

SLR1-8: 4 x 83mm holes

SLR2-14: 6 x 83mm holes

SLR2-8: 2 x 114mm holes

SLR1-22: 6 x 114mm holes



Polypropylene Spheres (Pack of 200)

Bath	NE1B-14	NE1B-18W	NE1B-28
No. of Packs	2 Packs		4 Packs
Accessory Code	BP0368		

NE1B Series Boiling Baths

Accessories



Raise and Lower Sterilising Tray

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	N/A	ST-1	N/A

Please note new models now come with ST-1 fitted as standard



Pair of guide rail assemblies for ST-1 tray

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	N/A	SA01409	N/A

Please note new models now come with SA01409 fitted as standard



Stainless Steel Test Tube Racks

Bath	NE1B-14	NE1B-18W	NE1B-28
No. of Racks	4		6
Accessory Code	6870 : 26 holes x 17mm DIA		
	6871 : 16 holes x 26mm DIA		
	6872 : 36 holes x 13mm DIA		
	6873 : 18 holes x 19mm DIA		
	6900 : 12 holes x 32mm DIA		



Drain Tap

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	N/A	BX0953	

Universal Accessories

CAL	Calibration Record
BX0616	Draining Syphon



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.