

# NE1B Analogue Boiling Bath Boiling Baths

*Clifton Range*<sup>®</sup>



NE1B-14

Designed for applications that require continuous boiling. This bath also has the flexibility where the control can be regulated from a gentle simmer to a vigorous boil.

Comprising of a stainless steel tank housed in a durable paint finished outer case. Painted surfaces feature a tough anti-bacterial finish with good chemical resistance. Complete with a removable stainless steel perforated false base. The bath, which is simmer-statically controlled, is ideally suited to prolonged periods of boiling.

Fitted with a sheathed immersion element. With the heater directly immersed in water it provides rapid boiling. A standard feature is a constant level device, which maintains a fixed level of liquid thereby eliminating the need to regularly top up the bath.




## NE1B - 14

## NE1B - 18W

## NE1B - 28

	NE1B - 14	NE1B - 18W	NE1B - 28
<b>GENERAL</b>			
Temp. Range	Simmer to boiling at 100°C		
Temp. Control	Analogue dial, energy regulator		
Heater Power (220-240V)	2500W		
Water Level Alarm	Features low liquid level float device, heating automatically suspended		
Voltage	230V / 50-60Hz		
Supplied with Lid		<input checked="" type="checkbox"/>	
Drain Tap Fitted		<input checked="" type="checkbox"/>	
<b>DIMENSIONS</b>			
Tank Capacity	14 Litre	18 Litre	28 Litre
External (mm)	362w x 332d x 290h		538w x 332d x 290h
Internal (mm)	325w x 300d x 150h	325w x 300d x 200h	500w x 300d x 200h

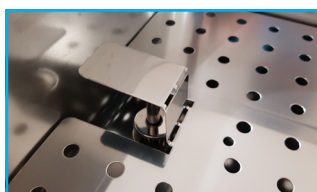

### Low liquid float device

A low liquid float device protects the heaters and a red warning light on the front panel alerts users to refill the bath.



### NE1B-18W

Raise and lower sterilising tray along with pair of guide rail assemblies comes as standard.



# NE1B Analogue Boiling Bath Accessories



## Stainless Steel Lift Off Gable Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	N/A	N/A	SL1-22SS



## Stainless Steel Flat Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	LD-14	N/A	LD-22



## Stainless Steel Hinged Gable Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	SL1-14H/SS	N/A	N/A



## Clear Hinged Gable Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	GL1-14		GL1-22



## Clear Lift Off Gable Lid (Fitted as standard)

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	GL1-14AF	N/A	GL1-22AF



## Stainless Steel Concentric Ringed Lid

Bath	NE1B-14	NE1B-18W	NE1B-28
Accessory Code	SLR1-14 SLR2-14	N/A	SLR1-22

SLR1-4: 2 x 114mm holes

SLR1-14: 4 x 114mm holes

SLR1-8: 4 x 83mm holes

SLR2-14: 6 x 83mm holes

SLR2-8: 2 x 114mm holes

SLR1-22: 6 x 114mm holes



## Polypropylene Spheres (Pack of 200)

Bath	NE1B-14	NE1B-18W	NE1B-28
No. of Packs	2 Packs		4 Packs
Accessory Code	BP0368		

# NE1B Analogue Boiling Bath Accessories



## Stainless Steel Test Tube Racks

Bath	NE1B-14	NE1B-18W	NE1B-28
No. of Racks	4		6
Accessory Code	6870 : 26 holes x 17mm DIA		
	6871 : 16 holes x 26mm DIA		
	6872 : 36 holes x 13mm DIA		
	6873 : 18 holes x 19mm DIA		
	6900 : 12 holes x 32mm DIA		

## Universal Accessories

CAL	Calibration Record
BX0616	Draining Syphon



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.