


Cole-Parmer® HM-200 and SHM-200 3 Position Multi-Mantles

- Individual built-in solid state electronic controls enable easy regulation of heater, while removing sparking associated with mechanical switching
- Three 12.7mm diameter support rods are included
- Round bottom flask capacity from 100mL to 1000mL at all 3 positions
- Top cover is polypropylene and external surfaces are powder coated giving good general chemical resistance
- Coiled heating element is suspended within a thermal insulating cartridge to provide maximum heat transfer and support
- Neon lights indicate when power is on and supplied to heaters and, if applicable, supplied to stirrers

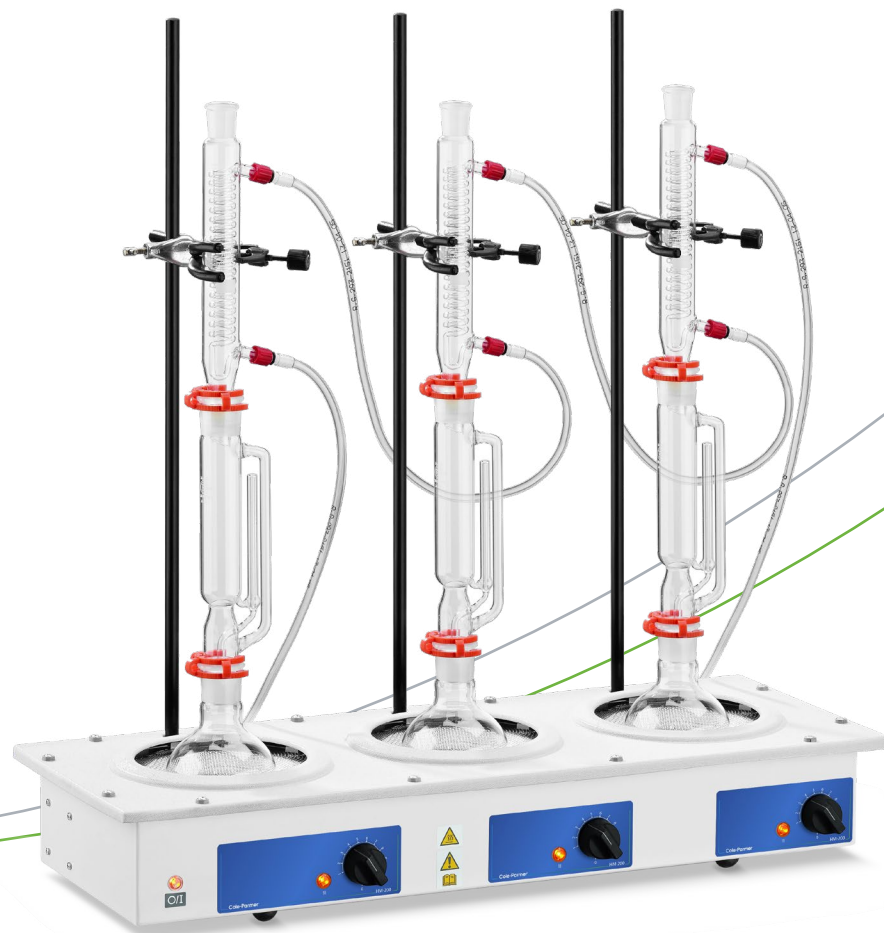

Ordering Information

Description	Ordering Number	Series No.	Model No.	Legacy Sku.
3 x 100ml, 230v 50/ 60Hz, 180W UK and EU plug	03076-01	HM-200	HM-200-MP3-100	EME30100/CEB
3 x 100ml, 115v 50/ 60Hz, 180W	03076-00	HM-200	HM-200-MP3-100-115	EME30100/CEBX1
3 x 250ml, 230v 50/ 60Hz, 450W UK and EU plug	03076-04	HM-200	HM-200-MP3-250	EME30250/CEB
3 x 250ml, 115v 50/ 60Hz/ 450W	03076-03	HM-200	HM-200-MP3-250-115	EME30250/CEBX1
3 x 500ml, 230v 50/ 60Hz, 600W UK and EU plug	03076-07	HM-200	HM-200-MP3-500	EME30500/CEB
3 x 500ml, 115v 50/ 60Hz, 600W	03076-06	HM-200	HM-200-MP3-500-115	EME30500/CEBX1
3 x 1000ml, 230v 50/ 60Hz, 900W UK and EU plug	03076-10	HM-200	HM-200-MP3-1000	EME31000/CEB
3 x 1000ml, 115v 50/ 60Hz, 900W	03076-09	HM-200	HM-200-MP3-1000-115	EME31000/CEBX1

3 x 100ml, 230v 50/ 60Hz, 220W UK and EU plug	03076-13	SHM-200	SHM-200-MP3-100	EMEA30100/CE
3 x 100ml, 115v 50/ 60Hz, 220W	03076-12	SHM-200	SHM-200-MP3-100-115	EMEA30100/CEX1
3 x 250ml, 230v 50/ 60Hz, 490W UK and EU plug	03076-16	SHM-200	SHM-200-MP3-250	EMEA30250/CE
3 x 250ml, 115v 50/ 60Hz, 490W	03076-15	SHM-200	SHM-200-MP3-250-115	EMEA30250/CEX1
3 x 500ml, 230v 50/ 60Hz, 640W UK and EU plug	03076-19	SHM-200	SHM-200-MP3-500	EMEA30500/CE
3 x 500ml, 115v 50/ 60Hz, 640W	03076-18	SHM-200	SHM-200-MP3-500-115	EMEA30500/CEX1
3 x 1000ml, 230v 50/ 60Hz, 940W UK & EU plug	03076-22	SHM-200	SHM-200-MP3-1000	EMEA31000/CE
3 x 1000ml, 115v 50/ 60Hz, 940W	03076-21	SHM-200	SHM-200-MP3-1000-115	EMEA31000/CEX1

Technical Specification

Specification	HM-200-MP3-100	HM-200-MP3-100-115	HM-200-MP3-250	HM-200-MP3-250-115	HM-200-MP3-500	HM-200-MP3-500-115	HM-200-MP3-1000	HM-200-MP3-1000-115
Flask Capacity	3 x 100ml	3 x 100ml	3 x 250ml	3 x 250ml	3 x 500ml	3 x 500ml	3 x 1000ml	3 x 1000ml
Dimensions (d x w x h), mm	270 x 640 x 110mm	270 x 640 x 110mm	270 x 640 x 110mm	270 x 640 x 110mm	285 x 755 x 130mm	285 x 755 x 130mm	285 x 755 x 130mm	285 x 755 x 130mm
No. of Recesses	3	3	3	3	3	3	3	3
Shipping weight, kg	6.2kg	6.2kg	6.2kg	6.2kg	7.4kg	7.4kg	7.4kg	7.4kg


Specification	SHM-200-MP3-100	SHM-200-MP3-100-115	SHM-200-MP3-250	SHM-200-MP3-250-115	SHM-200-MP3-500	SHM-200-MP3-500-115	SHM-200-MP3-1000	SHM-200-MP3-1000-115
Flask Capacity	3 x 100ml	3 x 100ml	3 x 250ml	3 x 250ml	3 x 500ml	3 x 500ml	3 x 1000ml	3 x 1000ml
Dimensions (d x w x h), mm	270 x 640 x 110mm	270 x 640 x 110mm	270 x 640 x 110mm	270 x 640 x 110mm	285 x 755 x 130mm	285 x 755 x 130mm	285 x 755 x 130mm	285 x 755 x 130mm
No. of Recesses	3	3	3	3	3	3	3	3
Shipping weight, kg	7kg	7kg	7kg	7kg	8.4kg	8.4kg	8.4kg	8kg

Cole-Parmer® HM-200 & SHM-200 3 Position Multi-Mantles

The Cole-Parmer® 3 place multi position (extraction) mantle has a built-in electronic controller and offers a choice of heater only (HM-200) or heater/stirrer (SHM-200) combination.

Additional Features for SHM-200 series models

- Allow the solution to be stirred and heated simultaneously
- Provide stirring speeds of between 50 to 1000rpm
- One stirring control operates each set of three recesses

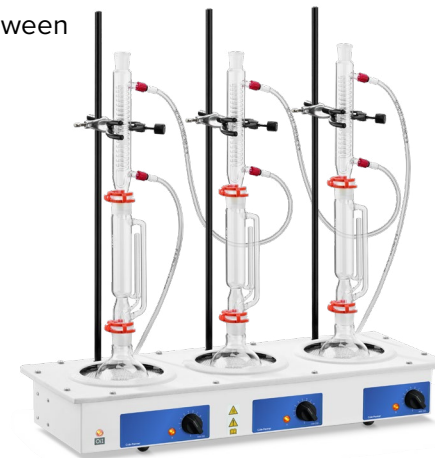



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.