

SK-R330-Pro

LCD Digital Rocking Shaker

SK-R1807-S

LED Digital Rocking Shaker

Features


- Rocking motion angle 9°(R330-Pro) and 7°(R1807-S)
- Two options:
10kg (R330-Pro) and speed range 10-70rpm,
2kg (R1807-S) and speed range 10-80rpm
- Double LCDs/LEDs for independent display of time and speed
- Continuous or timer operation mode
- Electronic time switching clock: 1min-99h59min
- Wide range of platforms for use with a variety of vessels(R330-Pro)
- Fixed work plate(R1807-S)
- Simple appearance, save operating space

Applications

Designed for time-dependent and low rpm sensitive shaking applications such as ELISA washes, enzyme immunoassays, protein synthesis, hybridization staining/destaining gels, immune precipitation and Western blots, cell culture, pharmaceutical profiling, etc.



SK-R330-Pro

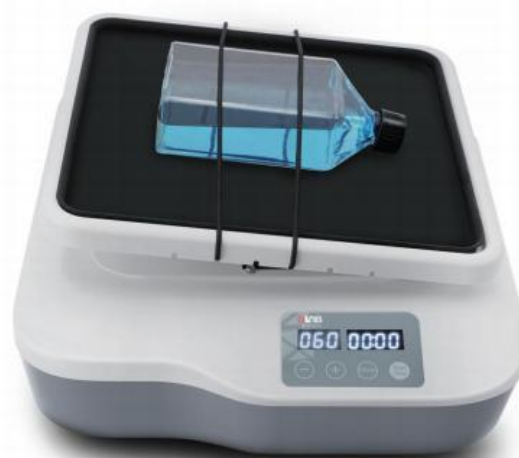

SK-R1807-S



Rocking Motion

Specifications	SK-R330-Pro	SK-R1807-S
Tilt angle	9°	7°
Max. load capacity(with platform)	10kg	2kg
Motor type	DC motor	DC motor
Motor input power	40W	20W
Motor output power	24W	16W
Speed range	10-70rpm	10-80rpm
Speed display	LCD	LED
Timer display	LCD	LED
Time settings range	1min-99h59min	1min-99h59min
Voltage,Frequency	100-240V, 50/60Hz	100-240V, 50/60Hz
Power	50W	20W
Weight	9kg	3.1kg
Permissible ambient temperature and humidity	5-40°C,80%RH	5-40°C,80%RH
Dimension[W×D×H]	360×410×200mm	294×348×123mm
Protection class	IP21	IP21

SK-D3309-Pro

SK-D1810-S

LCD Digital 3D Shaker

Features


- Provide a gentle three-dimensional motion shaker at a fixed tilt angle
- Anti slip silicone disc surface, better anti slip effect, and better sample processing effect
- Adjustable rotational mixing speed
- Independent LCD/LED display of time and speed
- Incubator compatible
- Various platform accessories (D3309-Pro)

Applications

3D shaker used for uniform shaking mixing in flasks, culture flasks, petri dishes and tubes. Ideal for culturing a variety of cell cultures, DNA extractions, distributing of cells and gel staining & destaining of cells.



SK-D3309-Pro

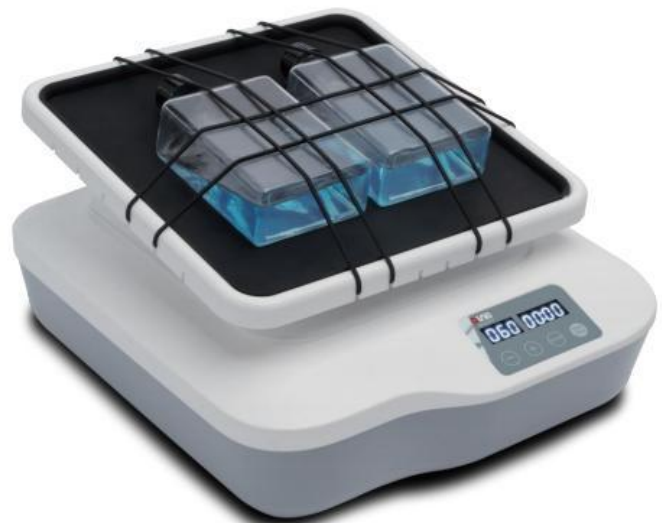

SK-D1810-S
















Wide range of accessories



3D Motion

Specifications	SK-D3309-Pro	SK-D1810-S
Tilt angle	9°	10°
Max. load capacity(with platform)	5kg	3kg
Motor type	DC motor	Brushless DC motor
Motor input power	40W	16W
Motor output power	24W	10W
Speed range	10-70rpm	10-80rpm
Speed display	LCD	LED
Timer display	LCD	LED
Time settings range	1min-99h59min	1min-99h59min
Voltage,Frequency	100-240V, 50/60Hz	100~240V,50/60Hz
Power	50W	20W
Weight	9kg	3.7kg
Permissible ambient temperature and humidity	5-40°C,80%RH	5-40°C,80%RH
Dimension[W×D×H]	360×430×106mm	294X348X170mm
Protection class	IP21	IP20



Cat. No.	Description	SK-O330-Pro	SK-L330-Pro	SK-R330-Pro	SK-D3309-Pro	SK-O330-M	SK-O180-Pro	SK-L180-Pro	SK-O180-S Multifunctional models
 18900155	SK330.5, Dish platform with non-slip mat and fixed string (8pcs) Dimension[WxD]: 30x28.5cm	■	■	■	■				
 18202835	Silicone mattress for 16× 5-15mL tubes, use with SK330.5	■	■	■	■				
 18202846	Silicone mattress for 6× 50mL tubes, use with SK330.5	■	■	■	■				
 18900185	SK180.4, Dish attachment with non-slip mat Dimension[WxD]: 24x21.5cm						■	■	■
 18900869	SK330.4 Dish platform with non-slip mat and fixed string(8 pcs), 34X34cm					■			
 18900859	SK330.6 Transparent shield, serologic test for syphilis, 41X39.5X10cm	■	■			■			
 18900860	SK180.5 Transparent Cover for Syphilis Serum Test, 32X32X8.3cm						■	■	
 18202852	Silicone mattress for 10× 5~10mL tubes						■	■	■
 18202853	Silicone mattress for 4× 50mL tubes						■	■	■
 18900027	SK330.1, Universal attachment with 4 bars Size: 32x32cm Capacity for flask: 25mLx21pcs, 50mLx15pcs, 100mLx12pcs, 250mLx6pcs, 500x4pcs	■	■	■	■	■			
 18900040	SK330.3 Lengthways roller attachment, Mainly used for separation funnel oscillation extraction applications 500mL*3pcs 250mL*4pcs 100/125mL*4pcs 50mL*5pcs	■	■	■	■	■			
 18900025	SK180.1, Universal attachment with 3 bars Size: 25x24.5cm Capacity for flask: 25mLx10pcs, 50mLx8pcs, 100mLx6pcs, 250mLx4pcs, 500x2pcs						■	■	■
 18900038	SK180.3 Lengthways roller attachment, Mainly used for separation funnel oscillation extraction applications 500mL*2pcs 250mL*3pcs 100/125mL*3pcs 50mL*4pcs						■	■	■

Cat. No.	Description	SK-O330-Pro	SK-L330-Pro	SK-R330-Pro	SK-D3309-Pro	SK-O330-M	SK-O180-Pro	SK-L180-Pro	SK-O180-S Multifunctional models
 18900036	Bar for SK330.1 universal attachment	■	■	■	■				
 18900037	Bar for SK180.1 universal attachment						■	■	■
 18900280	Connect sticks for double dishes shaker	■	■	■	■	■	■	■	■
 18900028	SK330.2, Fixing clip attachment Size: 32 x 32cm	■	■	■	■				
 18900026	SK180.2, Fixing clip attachment Size: 25 x 24cm						■	■	■
 18900029	Fixing clip for flask volume 25mL, used with SK330.2[9pcs] SK180.2[28pcs]	■	■	■	■		■	■	■
 18900030	Fixing clip for flask volume 50mL, used with SK330.2[9pcs] SK180.2[16pcs]	■	■	■	■		■	■	■
 18900031	Fixing clip for flask volume 100mL, used with SK330.2[9pcs] SK180.2 [16pcs]	■	■	■	■		■	■	■
 18900032	Fixing clip for flask volume 200/250mL, used with SK330.2[4pcs] SK180.2[9pcs]	■	■	■	■		■	■	■
 18900033	Fixing clip for flask volume 500mL, used with SK330.2[2pcs] SK180.2[9pcs]	■	■	■	■		■	■	■



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.