

## Genlab Multi-Purpose Benchtop Incubators



Model Shown - INC/100/TDIG/F



Model Shown - MINI/50/TDIG/F

The Genlab range of Benchtop Incubators offers a selection of highly efficient, reliable and cost effective incubators to suit most biological analysis, research and general laboratory applications. Available in a variety of sizes and designed with both side or base mounted displays to suit your exact needs. They feature our latest touch screen control system which offers intuitive control and excellent accuracies designed for laboratory incubators.

### Features

- 10 sizes 6 to 200 litres  
(base mount controller design)
- 6 sizes 50 to 200 litres  
(side mount controller design)
- Temp range (ambient+ 5) °C to 70 °C
- Fluctuation +/- 0.25°C @ 37 °C
- High resolution 0.1°C
- On screen historical trending (48hours)
- Manual or automatic overheat reset
- High accuracy Pt100B duplex sensors <0.8 °C
- Excellent stability <+/- 0.6 °C
- Low chamber uniformity

### Options

- Fan assisted circulation (not available <30L)
- Stainless steel chamber
- Viewing windows
- Inner glass doors
- Access ports (25, 50, 75 or 100mm)
- Bespoke stands and stacking kits
- Audible warnings
- 8 stage profile control
- Count down timer
- 24/7 real time clock operation
- Traceable calibration to national standards
- Extractor unit
- Extended warranty
- Internal lights
- Bespoke solutions available upon request



## Design


The exterior is constructed from sheet steel and finished in an easy clean powder coated paint. The interior chamber is made from mild steel coated with aluminium (CLAD), with a stainless steel chamber available as an option. They are fitted with fixed shelf runners and removable chrome plated wire grid shelves.

## Heating

Heated by Incoloy sheathed elements positioned in the base of the oven on all natural convection ranges. Fan assisted ranges have elements on the rear on base mount designs and in the side on side mount designs.

## Controls

The control system comprises of a bespoke touch screen user interface that offers both accurate temperature control and an integral overheat system by using 2 individual PT100 sensors. Advanced features can be turned on and off to suit bespoke heating applications and tailored to customer preference. Oven trending is displayed for up to 48 hours and optional upgrades are available for alarm outputs.

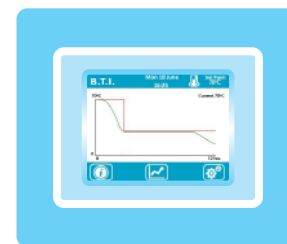

## Sizes and Specifications

Genlab Reference Number	Capacity* (litres)	Natural Convection		Fan Assisted Circulation		Shelves /Positions	Wattage	Weight (kg)
		Internal Dims. (H x W x D cms)	External Dims. (H x W x D cms)	Internal Dims. (H x W x D cms)	External Dims. (H x W x D cms)			
<b>Base Mount Designs...</b>								
MINI/6/TDIG	6	15 x 23 x 19	39 x 32 x 33	N/A	N/A	1..1	150	16
MINI/18/TDIG	18	26 x 26 x 27	52 x 38 x 41	N/A	N/A	1..1	350	21
MINI/30/TDIG	30	24 x 36 x 35	56 x 50 x 49	24 x 36 x 31	56 x 50 x 62	2..2	500	30
MINI/40/TDIG	40	32 x 36 x 35	64 x 50 x 49	32 x 36 x 31	64 x 50 x 62	2..3	500	35
MINI/50/TDIG	53	33 x 49 x 33	65 x 63 x 47	33 x 49 x 29	65 x 63 x 59	2..3	500	45
MINI/75/TDIG	73	33 x 49 x 45	65 x 63 x 60	33 x 49 x 41	65 x 63 x 72	2..3	750	55
MINI/100/TDIG	99	45 x 49 x 45	77 x 63 x 60	45 x 49 x 41	77 x 63 x 72	3..4	750	60
MINI/125/TDIG	121	55 x 49 x 45	87 x 63 x 60	55 x 49 x 41	87 x 63 x 72	3..5	750	65
MINI/150/TDIG	143	55 x 49 x 53	87 x 63 x 68	55 x 49 x 49	87 x 63 x 80	3..5	750	70
MINI/200/TDIG	180	N/A	N/A	75 x 49 x 49	107 x 63 x 80	4..7	1000	80
<b>Side Mount Designs...</b>								
INC/50/TDIG	58	42 x 40 x 35	58 x 73 x 48	45 x 36 x 35	58 x 73 x 48	2..3	500	35
INC/75/TDIG	78	42 x 53 x 35	58 x 85 x 48	45 x 49 x 35	58 x 85 x 48	2..3	750	45
INC/100/TDIG	102	42 x 53 x 46	58 x 85 x 59	45 x 49 x 40	58 x 85 x 59	3..4	750	55
INC/125/TDIG	126	52 x 53 x 46	69 x 85 x 59	55 x 49 x 46	69 x 85 x 59	3..5	750	65
INC/150/TDIG	149	52 x 53 x 54	69 x 85 x 67	55 x 49 x 54	69 x 85 x 67	3..5	750	70
INC/200/TDIG	198	N/A	N/A	75 x 49 x 54	89 x 85 x 67	4..7	1000	90

\* Note, Useable internal capacity is reduced with fan assisted circulation

## Ordering

All units have CLAD interior with touch screen interface as standard.

Where options are required simply add the option code to the original reference.

E.g. MINI/50/TDIG with fan assisted circulation and stainless steel chamber becomes MINI/50/TDIG/F/SS.

Contact sales for bespoke options.

## Option Codes

Fan Assisted Circulation - F

Double Door (Glass) - DD

Audible Warning - AWA

External Extraction System - EXT

Stainless Steel Chamber - SS

Access Ports - AP

3 Stage Advanced Timer - TIM

Extended Warranty - X1

Viewing Window - VIS

Stands/Stacking Kits - STACK

Traceable Calibration - SPC

Internal Lights - IL





# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.