



Multiparameter Waterproof Meters

Multiparameter measurement in compact,
IP67 rated waterproof enclosures

 **HANNA**[®]
instruments

HI98195

Multiparameter Waterproof Meter

pH, ORP, EC, TDS, Resistivity, Salinity, Seawater σ and Temperature

- **Waterproof**
 - IP67 rated waterproof, rugged enclosure for meter, IP68 for probe
- **Digital probe**
 - Digital probe with integral temperature sensor and two ports for pH (ORP) and EC sensors
- **Color coded, field replaceable sensors**
- **Auto-sensor recognition**
- **Quick calibration feature**
- **Automatic temperature compensation**
- **Automatic logging**
 - Store up to 45,000 samples
- **Log-on-demand**
 - Store measurement data at the press of a button
- **GLP**
 - GLP data provides data from previous five calibrations to ensure Good Laboratory Practices are met
- **Dedicated help key**
 - On-screen context specific help is readily available at the press of a button
- **Backlit LCD display with multifunction virtual keys**
- **Intuitive keypad**
 - Dedicated and virtual soft keys
- **Connectivity**
 - PC connectivity via opto-isolated micro-USB with Hanna software
- **Up to 360 hours of battery life**
 - Powered by (4) 1.5V AA batteries



- **Quick Connect Probe**
 - Fast, secure and waterproof connection

For Universal Applications

HI98195 provides multiparameter measurement in a compact and rugged, IP67 waterproof enclosure. Ideal for demanding applications, this meter features our rugged, quick connect multi-function probe with field replaceable sensors.

A backlit, graphic LCD provides easy to read resolution even in low-lit areas. A combination of dedicated and soft keys allows easy, intuitive operation in a choice of languages.

Data Logging

The HI98195 allows users to store up to 45,000 continuous or log-on-demand samples with logging intervals from one second to three hours. Logged data can be later transferred to a PC with the HI920015 micro USB cable and Hanna HI9298194 software.

GLP and On-Screen Help

Comprehensive GLP data includes data from last five calibrations. The contextual HELP menu can be accessed to obtain on-screen information and assistance about each feature at the touch of a button.


Quick Connect Probe

The HI7698195 probe features a quick connect DIN connector to make attaching and removing the probe simple and easy.




HI9828-25
Quick Calibration
solution

- **Field ready with Quick Calibration solution**
 - Standardize pH and conductivity with one calibration solution
 - Simply screw the calibration beaker filled with HI 9828-25 (500 mL) or HI9828-27 (1 gallon) solution onto the probe, select "Quick calibration" from the menu and press OK. Individual calibration may also be performed using multiple calibration points

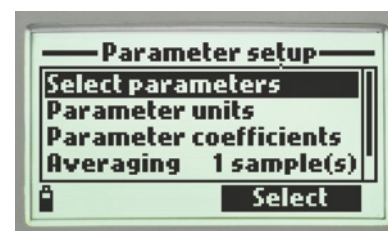

- **Setup**
 - Extensive setup screen features

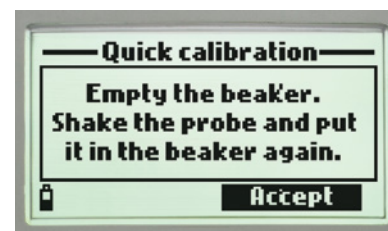

- **Guides**
 - On screen guides are displayed to help you while performing procedures such as calibration

pH: Specific Features

- **Calibration**
 - Up to a three-point calibration with five standard buffers and one custom buffer available
- **pH in mV option**
 - Useful for diagnostics
- **GLP data**
 - Offset, slope, date, time and buffers used
- **Automatically temperature compensated readings**
- **pH or pH/ORP field replaceable sensors**
 - Gel-filled and maintenance free
 - Double junction for reduced contamination of reference cell

EC/TDS/Resistivity: Specific Features

- **Calibration**
 - Single-point calibration from six standards (84 $\mu\text{S}/\text{cm}$, 1,413 $\mu\text{S}/\text{cm}$, 5.00 mS/cm , 12.88 mS/cm , 80 mS/cm , 111.8 mS/cm or custom point)
- **Temperature compensation**
 - Automatic Temperature Compensation
 - Configurable temperature coefficient range from 0.00 to 6.00%/°C
 - Choice of reference temperatures at 20 or 25°C
 - Absolute conductivity can be displayed along with the temperature compensated value
- **Auto-ranging**
- **Salinity readings**
 - Practical Salinity Scale (PSU) based on conductivity calibration



- **Rugged carrying case with custom insert included**

HI98195 Probe and Sensor Features and Specifications



| Specifications | HI7698195 | |
|----------------------------|--|--------------|
| Sensor Inputs | two (pH or pH/ORP, EC) | |
| Sample Environment | fresh, brackish, seawater | |
| Waterproof Protection | IP68 | |
| Operating Temperature | -5 to 55°C | |
| Storage Temperature | -20 to 70°C | |
| Maximum Depth | 20 m (66') | |
| Dimensions (without cable) | 342 mm (13.5"); 46 mm (1.8") dia | |
| Weight (without sensors) | 570 g (20.1 oz.) | |
| Cable Specification | multistrand-multiconductor shielded cable with internal strength member rated for 68 kg (150 lb.) intermittent use | |
| Wetted Materials | Body | ABS |
| | Threads | Nylon |
| | Shield | ABS / 316 SS |
| | Temperature Probe | 316 SS |
| | O-rings | EPDM |

Multi-function Sensor

- **Quick sensor replacement**
 - Sensor replacement is quick and easy with field replaceable, screw type connectors and are color coded for easy identification. These meters automatically recognize sensors



| Sensor Specifications | HI7698194-0 | HI7698194-1 | HI7698194-3 | |
|-----------------------|-----------------------------|---|---|-------------------------------------|
| Description | pH sensor | pH/ORP sensor | EC sensor | |
| Measurement Type | pH, mV (pH) | pH, mV (pH), ORP | EC | |
| Measurement Range | 0.00 to 13.00 pH; ±600.0 mV | 0.00 to 13.00 pH; ±600.0 mV; ±2000.0 mV | 0.0 to 200.0 mS/cm; 0.0 to 400 mS/cm (absolute) | |
| Temperature Range | -5 to 55°C | -5 to 55°C | -5 to 55°C | |
| Color Code | red | red | blue | |
| Materials | Tip | glass (pH) | glass (pH); Pt (ORP) | stainless steel electrodes AISI 316 |
| | Junction | ceramic | ceramic | |
| | Body | PEI | PEI | ABS/epoxy |
| | Electrolyte | gel | gel | |
| | Reference | double | double | |
| Maintenance Solution | HI70300 (storage solution) | HI70300 (storage solution) | none | |
| Dimensions | 118 x 15 mm | 118 x 15 mm | 111 x 17 mm | |
| Depth | 20 m (65') | 20 m (65') | 20 m (65') | |

| Specifications | HI98195 | |
|---------------------------|--|---|
| pH / mV | Range | 0.00 to 14.00 pH / ±600.0 mV |
| | Resolution | 0.01 pH / 0.1 mV |
| | Accuracy | ±0.02 pH / ±0.5 mV |
| | Calibration | automatic one, two, or three points with automatic recognition of five standard buffers (pH 4.01, 6.86, 7.01, 9.18, 10.01) or one custom buffer |
| ORP | Range | ±2000.0 mV |
| | Resolution | 0.1 mV |
| | Accuracy | ±1.0 mV |
| | Calibration | automatic at one custom point (relative mV) |
| EC | Range | 0 to 200 mS/cm (absolute EC up to 400 mS/cm) |
| | Resolution | manual: 1 µS/cm; 0.001 mS/cm; 0.01 mS/cm; 0.1 mS/cm; 1 mS/cm; automatic: 1 µS/cm from 0 to 9999 µS/cm; 0.01 mS/cm from 10.00 to 99.99 mS/cm; 0.1 mS/cm from 100.0 to 400.0 mS/cm; automatic mS/cm: 0.001 mS/cm from 0.000 to 9.999 mS/cm; 0.01 mS/cm from 10.00 to 99.99 mS/cm; 0.1 mS/cm from 100.0 to 400.0 mS/cm |
| | Accuracy | ±1% of reading or ±1 µS/cm whichever is greater |
| | Calibration | automatic single point, with six standard solutions (84 µS/cm, 1413 µS/cm, 5.00 mS/cm, 12.88 mS/cm, 80.0 mS/cm, 111.8 mS/cm) or custom point |
| TDS | Range | 0.0 to 400.0 ppt (g/L) (the maximum value depends on the TDS factor) |
| | Resolution | manual: 1 ppm (mg/L); 0.001 ppt (g/L); 0.01 ppt (g/L); 0.1 ppt (g/L); 1 ppt (g/L); automatic: 1 ppm (mg/L) from 0 to 9999 ppm (mg/L); 0.01 ppt (g/L) from 10.00 to 99.99 ppt (g/L); 0.1 ppt (g/L) from 100.0 to 400.0 ppt (g/L); automatic ppt (g/L): 0.001 ppt (g/L) from 0.000 to 9.999 ppt (g/L); 0.01 ppt (g/L) from 10.00 to 99.99 ppt (g/L); 0.1 ppt (g/L) from 100.0 to 400.0 ppt (g/L) |
| | Accuracy | ±1% of reading or ±1 ppm (mg/L) whichever is greater |
| | Calibration | based on conductivity or salinity calibration |
| Resistivity | Range | 0 to 999999 Ω•cm; 0 to 1000.0 kΩ•cm; 0 to 1.0000 MΩ•cm |
| | Resolution | dependent on resistivity reading |
| | Calibration | based on conductivity or salinity calibration |
| Salinity | Range | 0.00 to 70.00 PSU |
| | Resolution | 0.01 PSU |
| | Accuracy | ±2% of reading or ±0.01 PSU whichever is greater |
| | Calibration | based on conductivity calibration |
| Seawater σ | Range | 0.0 to 50.0 σ _v , σ ₀ , σ ₁₅ |
| | Resolution | 0.1 σ _v , σ ₀ , σ ₁₅ |
| | Accuracy | ±1 σ _v , σ ₀ , σ ₁₅ |
| | Calibration | based on conductivity or salinity calibration |
| Temperature | Range | -5.00 to 55.00°C; 23.00 to 131.00°F; 268.15 to 328.15K |
| | Resolution | 0.01°C; 0.01°F; 0.01K |
| | Accuracy | ±0.15°C; ±0.27°F; ±0.15K |
| | Calibration | automatic at one custom point |
| Additional Specifications | Temperature Compensation | automatic from -5 to 55°C (23 to 131°F) |
| | Logging Memory | 45,000 records (continuous logging or log-on-demand of all parameters) |
| | Logging Interval | one second to three hours |
| | PC Connectivity | via USB (with Hanna PC software) |
| | Environment | 0 to 50°C (32 to 122°F); RH 100% IP67 |
| | Battery Type / Life | 1.5V AA batteries (4) / approximately 360 hours of continuous use without backlight (50 hours with backlight) |
| Dimensions / Weight | 185 x 93 x 35.2 mm (7.3 x 3.6 x 1.4") / 400 g (14.2 oz.) | |



HI98194 includes:

HI98195 includes:



HI7698194/4
probe



HI7698194-0
pH sensor



HI7698194-3
EC sensor



HI7698194-2
DO sensor



HI7698195/4
probe



HI7698194-0
pH sensor



HI7698194-3
EC sensor



HI9828-25
quick calibration
solution



HI76981942
probe maintenance
kit



HI9828-25
quick calibration
solution



HI76981952
probe maintenance
kit



HI98196 includes:



HI7698196/4 probe



HI7698194-0 pH sensor



HI7698194-2 DO sensor



HI9828-25 quick calibration solution



HI76981942 probe maintenance kit



rugged carrying case with custom insert



HI76981943 calibration beaker



HI9298194 PC software



HI920015 micro USB cable



quality certificate



instruction manual



1.5V AA batteries (4)

All meters are also supplied with:



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.