



2xOF302-40MQ2 Oil Free Compressor with Dryer

JUN-AIR®



Features

- 100% oil free
- 100% duty cycle
- Low noise and vibration level
- Easy maintenance
- Long-life, 2-year warranty
- 8000 hours cup life at 8bar
- 4000 hours maintenance interval for the dryer

Description

- Built-in 2-column adsorption dryer
- Pressure dew point down to -40°C
- 1µm air filtration (class 2 for dust according to ISO8573.1:2001)
- Pressure regulator included
- Receiver with internal coating and 5-year warranty


Other Versions

- Without dryer: 2xOF302-40M
- Other power supply: 100V, 120V, 60Hz...

Product Specifications


Model	2xOF302-40MQ2		
Voltage	V	1x230	1x120
Frequency	Hz	50	60
Motor	HP	1.2	1.2
	kW	0.9	0.9
Displacement	l/min	216	276
Air flow @ 8bar	l/min	61	70
	cfm	2.15	2.5
Max. pressure	bar	8	8
	psi	120	120
Max. current	A	6.8	13.2
	litres	40	40
Receiver volume	gallon	10.6	10.6
	kg	112	112
Weight	lbs	247	247
	Noise level	dB(A)/1m	50
Duty cycle		100%	100%

Air Flow / Pressure Curve*

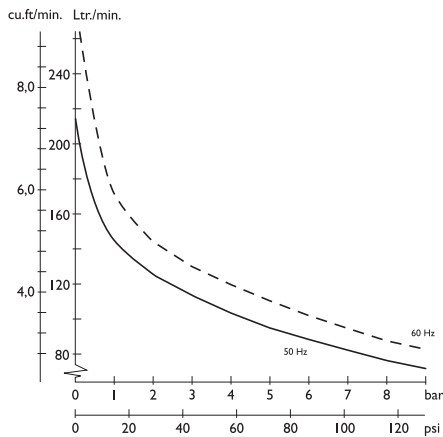

*air flow is reduced by approx. 18-20% on units with dryer

Dimensions*

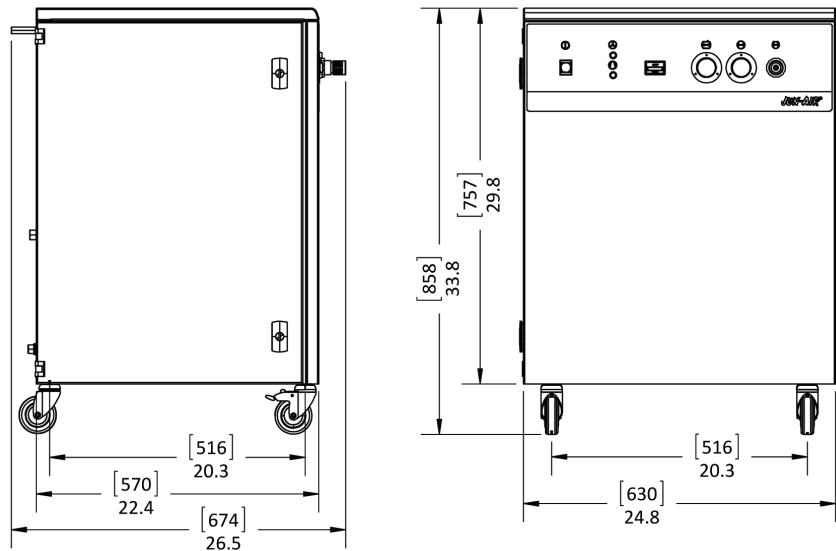

*for unit without accessories or options



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.