



# OF302-4MQ2 Oil Free Compressor with Dryer

# JUN-AIR®



### Features

- 100% oil free
- 100% duty cycle
- Low noise and vibration level
- Easy maintenance
- Long-life, 2-year warranty
- 8000 hours cup life at 8bar
- 4000 hours maintenance interval for the dryer

### Description

- Built-in 2-column adsorption dryer
- Pressure dew point down to -40°C
- 1µm air filtration (class 2 for dust according to ISO8573.1:2001)
- Pressure regulator included
- Receiver with internal coating and 5-year warranty


### Other Versions

- Without dryer: OF302-4M
- Other power supply: 100V, 120V, 60Hz...

### Product Specifications


Model	OF302-4MQ2		
Voltage	V	1x230	1x120
Frequency	Hz	50	60
Motor	HP	0.6	0.6
	kW	0.44	0.44
Displacement	l/min	108	138
Air flow @ 8bar	l/min	30	36
	cfm	1.06	1.27
Max. pressure	bar	8	8
	psi	120	120
Max. current	A	3.4	6.6
Receiver volume	litres	4	4
	gallon	1.1	1.1
Weight	kg	52	52
	lbs	115	115
Noise level	dB(A)/1m	47	48
Duty cycle		100%	100%

### Air Flow / Pressure Curve\*

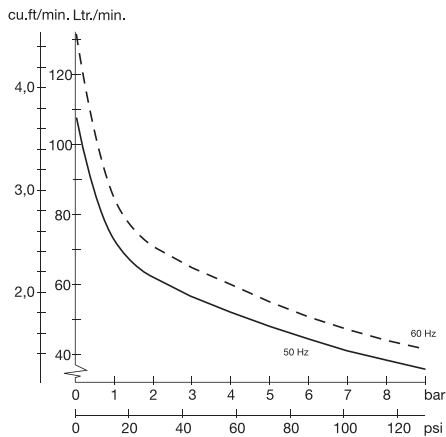

\*air flow is reduced by approx. 18-20% on units with dryer

### Dimensions\*

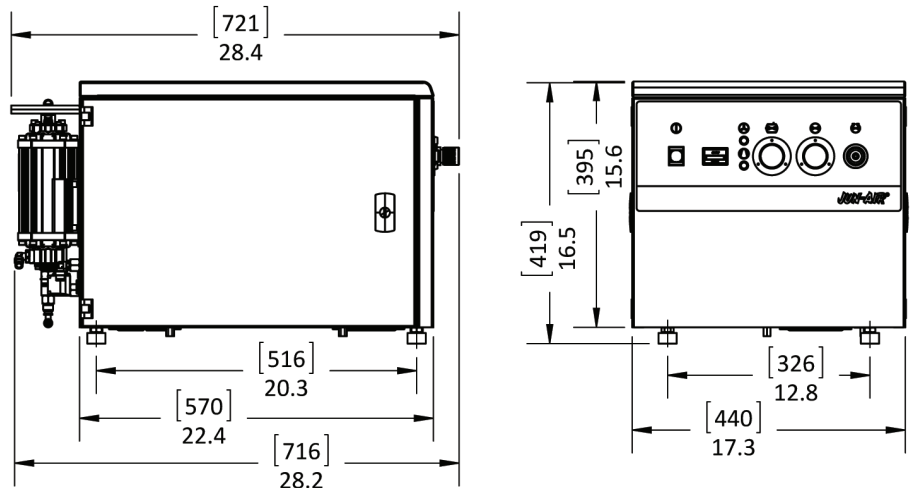

\*for unit without accessories or options



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.