


# 1 Door, Manual Close, Mini Transportable Flammable Cabinet for Aerosols, Sure-Grip® EX, Yellow - 890200

## Product Images



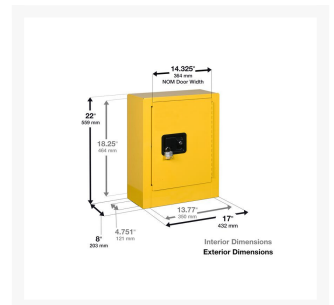

## Short Description

---

- Ideal when space is limited
- Protect workers
- Reduce fire risks
- Meet safety codes
- Ten-year warranty

## Description

---

Mini Safety Cabinets offer the smallest footprint for easy mobility on service vehicles, carts and maintenance trikes. Includes 3 bungee cords to keep contents in place during movement. These cabinets have all the great features and benefits of our classic cabinets as outlined below.

Protect workers, reduce fire risks, and improve productivity by storing flammable fuels and chemicals in code-compliant safety cabinets. Designed to meet OSHA and NFPA standards, Sure-Grip® EX cabinets are constructed of sturdy 18-gauge (1-mm) thick double-wall, welded steel with 1-1/2" (38-mm) of insulating air space for fire resistance. High-performance, manual-close, self-latching doors open to a full 180 degrees and close easily and securely for maximum protection under fire conditions when pushed shut. Fail-safe closing mechanism ensures three-point stainless steel bullet latching system works every time.

For added security, U-Loc™ handle comes with a cylinder lock and key set or accepts an optional padlock. Haz-Alert™ reflective warning labels are highly visible under fire conditions or during power outages when

illuminated with a flashlight.

Shelve adjust on 3" (76-mm) centers for versatile storage. Cabinets also feature a leakproof 2" (51-mm) bottom sump, four adjustable self-leveling feet, grounding connector and trilingual warning label. Durable lead-free epoxy/polyester powder-coat finish provides excellent chemical resistance.

All cabinets come with a ten-year warranty.

## Specifications

---

Model No	890200
UPC	697841145030
International Model No.	8902001
Color	Yellow
Pallet Quantity	8
Material Specifications	18-Gauge Steel
OSHA Compliance	Yes
Complies with NFPA Code 30	Yes
Number of Doors	1
Door Type	Manual
Number of Shelves	1
Adjustable Shelves	Yes
Dimensions, Exterior	22" H x 17" W x 8" D
	559mm H x 432mm W x 203mm D
Dimensions, Interior	18.25" H x 13.8" W x 4.75" D
	464mm H x 350mm W x 121mm D
Height (Inches)	22
Width (Inches)	17
Depth	8
Net Weight, lbs	23.0
Net Weight, kg	10.4



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.