

These multi-functional modular table systems are ideal for use within the medical, research and educational establishments.

Our Fusion Table Systems are fully compatible with our Fusion Unit Range allowing the user to easily change the functionality of their modular system in accordance with their own preference and needs.



Freestanding Table

These fully welded table frames come complete with a Trespa Athlon worktop and incorporates either adjustable feet or industrial castors.

The worktops have chamfered edges and curved corners. These tables are 1500 x 750 x 900mm high and support a load of 210kg spread evenly over the worktop. All steel framing is manufactured from 40x40x1.6mm

thick electric resistance welded seamless tubular steel box with all joints fully welded. Frames are then thoroughly de-greased prior to being Epoxy Powder Coated.



Mobile Table



Mobile Table with 3 No units

Each table is able to support up to 3 storage units. Units can be purchased at the same time as the tables or they can be easily installed into an existing bench from the Fusion Range at a later date.



Mobile Table with 2 No units

All these units are compatible with all the tables on the Fusion range and are designed to be suspended under the tables. Our unique suspension system allows the units to slide within the frame or be easily removed/added.

The overall sizes of the carcasses are 700mm High x 450mm Wide x 500mm Deep. Each unit weighs approximately 50kg, this has to be taken into account when loading the table.



Suspended 4 Drawer

Specification

The under bench storage units are manufactured from Melamine Faced MDF. The carcasses are manufactured using mitred joints and fully grooved in backs making for a long lasting and durable unit. All seen edges sealed with 1mm thick ABS lipping. The drawers and doors are also manufactured from Melamine faced MDF. All seen edges are sealed with 2mm thick ABS lipping's and come complete with a stainless steel 'D' Pull style handles.



Suspended Drawer & Cupboard



Suspended Cupboard

Kinetic

flexible solutions

Fusion Range of Modular tables and Suspended units 2014 price list

These multi-functional modular table systems are ideal for use within the medical, research and educational establishments.

Our Fusion Table Systems are fully compatible with our Fusion Unit Range allowing the user to easily change the functionality of their modular system in accordance with their own preference and needs.



Mobile Table



Freestanding Table




Suspended 4 Drawer



**Suspended Drawer
& Cupboard**



Suspended Cupboard

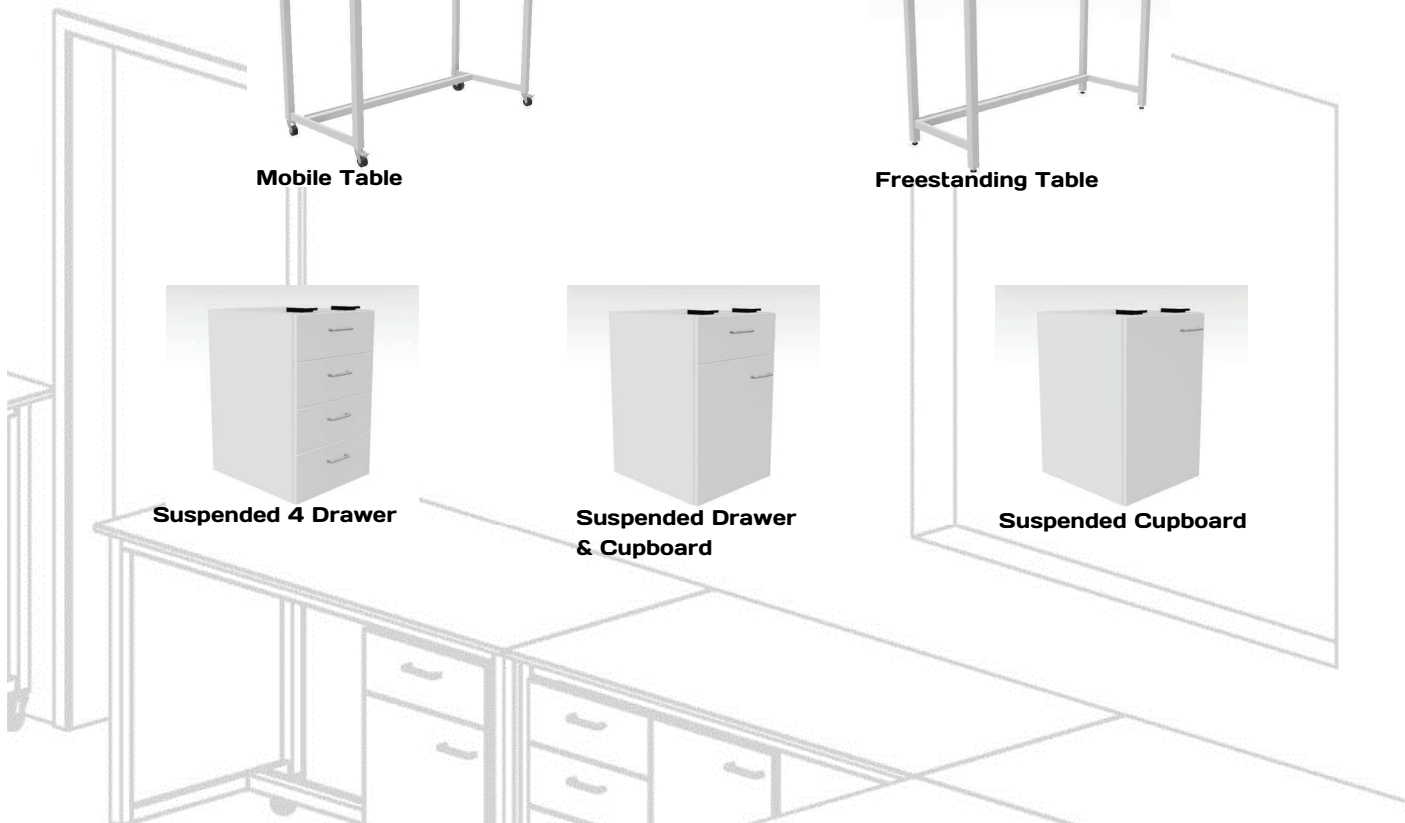



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.