




USB TEMPERATURE DATA LOGGER



- ✓ **Temperature Range of -35°C to +80°C.**
- ✓ **Simple and Reliable way of recording**
- ✓ **User can select recording intervals**
Every 10 seconds up to every 12 hours.
- ✓ **Can store up to 16, 382 readings.**
Which gives a memory of up to 22 years if set at the maximum log rate, and 3.5 months at the minimum.
- ✓ **Easy to download the results**
Simply plug the data logger into any USB port on your computer and download the stored data.
- ✓ **IP67 moisture resistant when cap is on**
So it is fine even if there is condensation, or a small spillage in the fridge or freezer.
- ✓ **User can select recording intervals**
Every 10 seconds up to every 12 hours.
- ✓ **Replaceable lithium battery**
- ✓ **Supplied with software CD, compatible with Windows 2000, XP, Vista and 7**

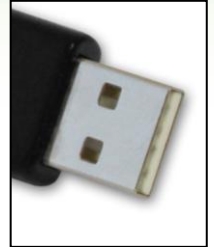

www.facebook.com/Lec-Medical
www.twitter.com/Lec_Medical

SKU CODE:

PRODUCT CODE:

DESCRIPTION:

444448729

ATMDL01

USB Temperature Data Logger




FULL SPECIFICATION – ATMDL01 Temperature Data Logger



Description	Data Logger
SKU code	444448729
GDPA code	ATMDL01
Unit External Width (mm)	25
Unit External Depth (mm)	25
Unit External Height (mm)	100
Packed Width (mm)	45
Packed Depth (mm)	105
Packed Height (mm)	115
Gross Capacity (L)	N/A
Net Capacity (L)	N/A
Type of plastic	N/A
Door Finish	N/A
Temperature Range	-35°C to +80°C
Temperature Accuracy	±1°C
IP67 moisture resistant	Yes (only when cap is on)
Supplied with software	Yes

EL-WIN-USB (Control Software) Features:

- Logger Name.
- °C or °F.
- Logging Rate (10s, 1m, 5m, 30m, 1h, 6h, 12h).
- High & Low Alarms
- Immediate, delayed & push-to-start settings.
- Display off, on for 30 seconds after button press, or permanently on.
- Data Rollover – *allows unlimited logging periods by overwriting the oldest data when the memory is full.*

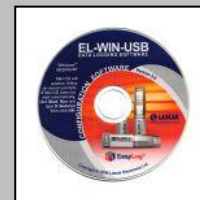

*The software compatible with Windows XP, Vista and 7 (including 64-bit operating systems). The software is **NOT** compatible with Mac or Linux.*



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.