

# Bluetooth Enabled Pharmacy Plus Range

47L Bluetooth Pharmacy Countertop Fridge Glass Door

GDPA code      PPGR47BT-UK  
SKU code        444411306




## Key features

- ✓ Local Monitoring and Alarm History via Bluetooth
- ✓ Power Failure Alarm with battery back up. Rechargeable battery
- ✓ Solid and glass door options
- ✓ Intelligent Fan Management for Enhanced Performance



### New Controller

Offers Local temperature monitoring and alarm history via Bluetooth. Power failure alarm with battery backup, rechargeable battery and intelligent fan management for enhanced performance

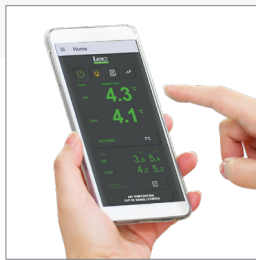

### New Bluetooth Enabled App

Accessibility features include contactless local monitoring of air and load temperatures. Downloadable temperature history and alarm history.



### Flexible internal storage

Multiple shelf locations along with a wire basket provides a flexible storage solution.



### Anti-microbial handle

Powder coated to help prevent cross contamination by inhibiting the growth of microbes.



### Dual probe

Greater temperature accuracy given from two temperature probe readings. One probe measuring air temperature and one measuring load (or product) temperature.

Specification: Upright 41L Bluetooth Pharmacy Countertop Fridge Glass Door

<b>Gross Capacity (litres)</b>	47L
<b>Net Capacity (litres)</b>	30L
<b>Unit External Width (mm)</b>	475mm
<b>Unit External Depth (mm)</b>	445mm
<b>Unit External Height (mm)</b>	565mm
<b>Product Weight (kg)</b>	17.5kg
<b>Packed Width (mm)</b>	490mm
<b>Packed Depth (mm)</b>	530mm
<b>Packed Height (mm)</b>	620mm
<b>Packed Weight (kg)</b>	21.5kg
<b>Country of Origin</b>	UK
<b>Internal Dimensions (mm) (WxDxH)</b>	403mm W x 354mm D x 425mm H
<b>Temperature Range</b>	2°C to 8°C
<b>Colour</b>	White/green
<b>External Finish</b>	Painted Steel
<b>Internal Finish</b>	HIPS*
<b>Type of plastic</b>	SRL800 (HIPS/ABS mix)
<b>Door Finish</b>	PVC/Glass
<b>Door Type</b>	Glass
<b>Door Style</b>	Hinged, reversible
<b>Door Liner</b>	Flat
<b>Door Open Depth (mm)</b>	905mm
<b>Replacement Door Seal</b>	No
<b>Open Door Sensor</b>	Yes
<b>Lock</b>	Yes, 2 Keys
<b>Adjustable Feet (number &amp; position)</b>	2 adjustable (front)
<b>Porthole</b>	Yes
<b>Porthole Size</b>	16 mm
<b>Porthole Location</b>	Rear LH Corner
<b>Temperature Display on Panel</b>	Yes
<b>Temperature Alarm</b>	Yes
<b>Remote Alarm Contacts</b>	Yes
<b>Audible Alarms</b>	Yes
<b>Visual Alarms</b>	Yes
<b>Number of Shelves</b>	1
<b>Number of Shelf Positions</b>	3
<b>Shelf Size (mm) (WxDxH)</b>	416mm W x 200mm D x 6mm H
<b>Shelf Material</b>	Polythene Coated Wire
<b>Number of Baskets</b>	1
<b>Basket Material</b>	Polythene Coated Wire
<b>Basket Size (mm) (WxDxH)</b>	396mm W x 165mm D x 150mm H
<b>Interior Light</b>	Yes
<b>Light Type</b>	LED
<b>Light Activation</b>	Door Swith / Controller
<b>Fan</b>	Yes - Axial
<b>Spark Free</b>	No
<b>Defrost setting</b>	Auto (Off Cycle)
<b>Controller Type</b>	Digital BT, external
<b>Number of Probes</b>	2
<b>Probe Type</b>	Air & Load (Glycol Solution)
<b>Ambient Temperature Range</b>	10°C to 32°C
<b>Refrigerant Type</b>	R600a
<b>CFC Free</b>	Yes
<b>HFC Free</b>	Yes
<b>HCFC Free</b>	Yes
<b>Ammonia Free</b>	Yes
<b>UK 3 Pin Plug</b>	Yes
<b>Other Plug Types Available</b>	Yes
<b>Voltage Supply</b>	220 -240V
<b>Insulation Thickness (mm)</b>	40mm
<b>Insulation Type</b>	High Density Polyurethane Foam
<b>Blowing Agent</b>	Cyclopentane
<b>Anti Microbial Handle</b>	Yes
<b>Wall Brackets</b>	Yes
<b>Optional Drawer System</b>	No
<b>Optional Digital Lock</b>	No
<b>Other Door Option</b>	Solid PPSR47UK - 444410570
<b>Battery Backup</b>	Yes



# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.