



Designed for Routine Weighing Applications in Your Workplace

Offering accuracy and repeatability in essential weighing applications in laboratory, industrial and education settings, PR Analytical and Precision Balances deliver competitive performance at an economical price. Featuring RS232 connectivity for easy communication, and a backlit display and a simple interface for uncomplicated operation, the PR is perfectly designed for your workplace.

Standard Features Include:

- **Basic Functionality for Routine Weighing Applications**

The PR is equipped with three essential weighing modes, RS232 connectivity for data transfer and storage and internal calibration (in select models), making it ideal for routine weighing applications.

- **Designed for Uncomplicated Operation with Easy-to-Use Display and Interface**

Equipped with an easy-to-read, bright backlit display and a simple user interface, the PR is incredibly easy to operate, with almost no training required.

- **Smart Design and Durable Construction**

The PR's small footprint saves desktop space, while providing a large weighing surface. The PR is durably constructed, and features a stainless steel pan to withstand day-to-day use in the workplace.

PR SERIES *Analytical and Precision Balances*

Standard RS232 connectivity for easy data transfer.



A large, bright backlit display allows for easy viewing of results.



Durably constructed, and features a stainless steel pan to withstand day-to-day use.

Two adjustable feet that allow user to level the balance.

Four mechanical keys enable easy operation of basic application modes.

InCal™ Model	PR124	PR224	PR223	PR423	PR523
ExCal Model	PR124/E	PR224/E	PR223/E	PR423/E	PR523/E
Approved Model*	PR124M	PR224M	PR223M	PR423M	PR523M
Capacity (g)	120	220		420	520
Readability (g)	0.0001		0.001		
Verification Interval (e) (g) (Approved Models)	0.001		0.01		
Class (Approved Models)	I		II		
Repeatability (STDEV) (g)	0.0001		0.001		
Linearity (g)	0.0002		0.002		
Stabilization Time (s)	3		2		
Sensitivity Temperature Drift (PPM/K)	±3		±8	±3	
Typical Minimum Weight USP (USP, K=2, U=0.10%)	200 mg		2 g		
Optimized Minimum Weight (USP, U=0.10%, K=2) SRP≤0.41d**	82 mg		0.82 g		
Units	Milligram, Gram, Ounce, Carat, Pennyweight, Ounce Troy, Newton, grain				
Units (Approved Models)	Milligram, Gram, Carat				
Applications	Basic Weighing, Parts Counting, Percent Weighing				
Pan Size (Ø)	90 mm		120 mm		
Tare Range	To full range by subtraction				
Power Supply	Power Input: 100 – 240V ~ 200mA 50 – 60Hz 12 – 18VA Power Output: 12 VDC 0.5A				
Assembled Dimensions (W × D × H)	201 × 317 × 303 mm				
Communication	RS232				
Operating Temperature Range	Operating conditions for ordinary lab application: +10 to 30°C (operability guaranteed between +5 and 40°C)				
Storage Temperature Range	Humidity: maximum relative humidity 80% for temperatures up to 30 °C, decreasing linearly to 50% relative humidity at 40 °C				
Net Weight	4.5 kg				
Shipping Weight	7 kg				
Shipping Dimensions (W × D × H)	507 × 387 × 531 mm				

*Approved models are all internal calibration models.

**The value for SRP is the standard deviation for n replicate weighings (n≥10)

PR SERIES *Precision Balances*

Standard RS232 connectivity for easy data transfer.



A large, bright backlit display allows for easy viewing of results.



Two adjustable feet that allow user to level the balance.

Four mechanical keys enable easy operation of basic application modes.




InCal™ Model	PR1602	PR2202	PR4202	PR4201	PR6201
ExCal Model	PR1602/E	PR2202/E	PR4202/E	PR4201/E	PR6201/E
Approved Model*	PR1602M	PR2202M	PR4202M	PR4201M	PR6201M
Capacity (g)	1600	2200	4200		6200
Readability (g)	0.01			0.1	
Verification Interval (e) (g) (Approved Models)	0.1				1
Class (Approved Models)	II				
Repeatability (STDEV) (g)	0.01			0.1	
Linearity (g)	0.02			0.2	
Stabilization Time (s)	1				
Sensitivity Temperature Drift (PPM/K)	±6		±3	±10	
Typical Minimum Weight USP (USP, K=2, U=0.10%)	20 g			200 g	
Optimized Minimum Weight (USP, U=0.10%, K=2) SRP≤0.41d**	8.2 g			82 g	
Units	Gram, Kilogram, Ounce, Pound, Carat, Pennyweight, Ounce Troy, Newton, grain				
Units (Approved Models)	Gram, Kilogram, Carat				
Applications	Basic Weighing, Parts Counting, Percent Weighing				
Pan Size (Ø)	180 mm				
Tare Range	To full range by subtraction				
Power Supply	Power Input: 100 – 240V ~ 200mA 50 – 60Hz 12 – 18VA Power Output: 12 VDC 0.5A				
Assembled Dimensions (W × D × H)	201 × 317 × 93 mm				
Communication	RS232				
Operating Temperature Range	Operating conditions for ordinary lab application: +10 to 30°C (operability guaranteed between +5 and 40°C)				
Storage Temperature Range	Humidity: maximum relative humidity 80% for temperatures up to 30 °C, decreasing linearly to 50% relative humidity at 40 °C				
Net Weight	3.5 kg				
Shipping Weight	5 kg				
Shipping Dimensions (W × D × H)	550 × 385 × 291 mm				

*Approved models are all internal calibration models.

**The value for SRP is the standard deviation for n replicate weighings (n≥10)

PR SERIES Analytical and Precision Balances

Dimensions

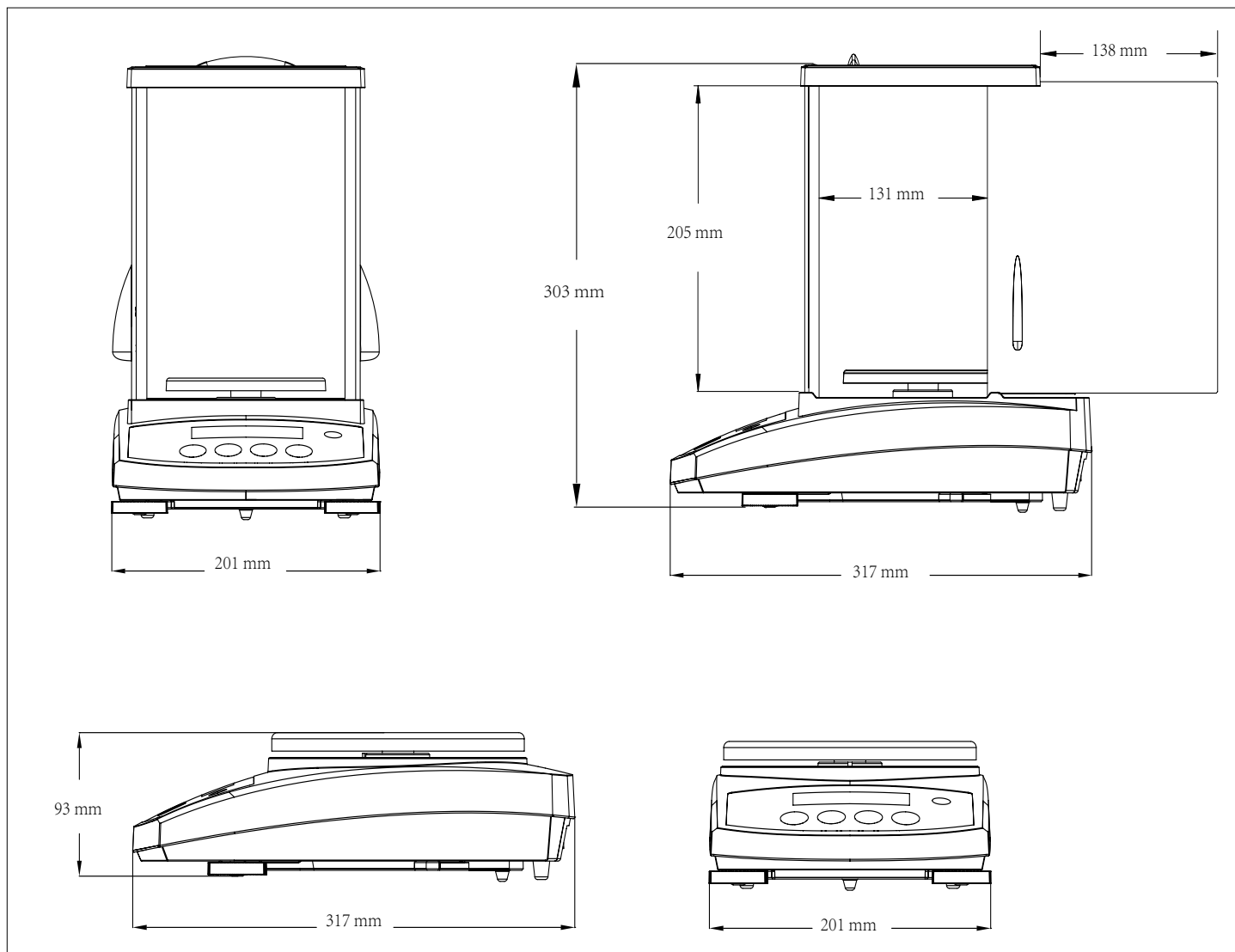

Other Standard Features and Equipment

ABS top housing, removable stainless steel pan, glass draftshield with sliding top door, integrated weigh-below-hook, security bracket, calibration lock, user-selectable environmental filters and brightness settings, auto-tare, auto-dim, user-selectable span calibration points, software lockout and reset menu, user-selectable communication settings and data print options, user-definable project and user IDs, software overload/underload indicator, stability indicator

Compliance

- **Metrology:** OIML R76; EN 45501
- **Product Safety:** IEC/EN 61010-1; CAN/CSA C22.2 61010-1; UL 61010-1
- **Electromagnetic Compatibility:** IEC/EN 61326-1 Class B, basic environments; FCC Part 15 Class A; Canada ICES-003 Class A
- **Compliance Marks:** CE; CSA

Accessories

Auxiliary Display AD7-RS	30472064
Density Kit for Solids	80253384
Sinker Glass for Density Determination	83034024
Printer SF40A	30064202
SF40A Ink Ribbon Cassette	12120798
SF40A Paper Roll 57.5mm (2 Rolls)	12120799
Security Device	80850043
In-Use Cover	30372547
Dust Cover (for 0,1mg and 1mg models)	30093334
Stand-alone Ionizer ION-100A	30095929
RS232 Cable (9-pin)	00410024

8077257_B 20180416 © Copyright OHAUS Corporation





WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.