

PURITE

Purite Neptune Ultimate

Reliable water purification system, consistently delivering ultrapure water for laboratories



The Purite Neptune Ultimate provides a high flow of guaranteed 18.2MΩcm ultrapure water for analytical and life science laboratory applications.

The Purite Neptune Ultimate utilises a number of proven technologies to produce ultra-pure, 18.2MΩ.cm water on-demand

Additional features

- Ion exchange cartridges incorporating monospherical, semiconductor grade mixed bed deionising resin with a low TOC leaching profile and a high activity organic absorption media
- Sub-micron and ultrafilter to remove particles, bacteria and endotoxins
- 185nm or photo-oxidising UV to cleave organic compounds into smaller charged ionic species that can be removed by ion exchange
- 254nm UV to reduce bacteria by more than 99%
- UV LED disinfection incorporated within the dispense head
- WiFi enabled for remote monitoring and operation.

Purite Neptune Ultimate


Reliable water purification system, consistently delivering ultrapure water for laboratories

Common Purite Neptune Ultimate features

Our Purite Neptune Ultimate is compact, robust, simple to use and easy to maintain. Other features of the Purite Neptune Ultimate include:

- Space-saving, bench top or wall mounted systems
- Dispense rate up to 2l/min
- Installation kit and all consumables included for first few months operation
- LCD colour touch screen panel
- Visual and audible alarms included
- USB port to download event data and upload software dates
- Semi-automatic clean cycle
- TOC monitoring.

Purite Neptune Ultimate Process Flow

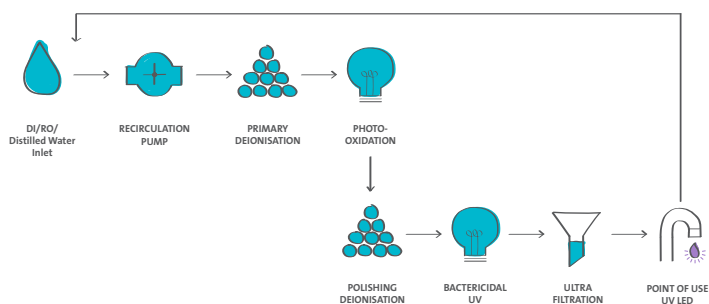

Technical Specifications

Unit Specification	
Width (mm)	310
Depth (mm)	560
Height (mm)	750
Max shipping weight (kg)	21
Max working weight (kg)	29
Installation requirements	
Power	Single Phase, 110-230V, +/- 10%, 50/60 Hz
Feed water	< 20µS/cm
Maximum TDS (ppm)	< 14
Minimum inlet pressure - psi (bar)	0.5 (Flooded Suction)
Maximum inlet pressure - psi (bar)	20 (1.38)
Feed water temperature	1-35°C

System Specification	
Display panel	LCD – Colour touch screen
Internal filtration	Ultrafiltration
Point of use	UV / LED / 0.2µm (Optional)
UV lamp	185nm & 254nm
Recirculation pump	✓
Ultrapure polishing cartridges	✓
Dispense modes	Latched - hold - volumetric
Dispense flow rate	up to 2.0 l/min

Treated Water Specification	
Inorganics	18.2MΩ.cm
pH ¹	Neutral
Bacteria	< 0.1cfu/ml
Organics – TOC (ppb)	< 1
Particles	0.005µm
Endotoxins	< 0.001EU/ml
DNases	< 4pg/µl
RNases	< 0.01ng/ml

* pH of stored water may decrease due to absorption of free carbon dioxide

PURITE



WolfLabs

Pricing on any accessories shown can be found by keying the part number into the search box on our website.

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

www.wolflabs.co.uk

Tel : 01759 301142

Fax : 01759 301143

sales@wolflabs.co.uk

Please contact us if this literature doesn't answer all your questions.