


# Waysafe SA640

## Sample Analysis Enclosure

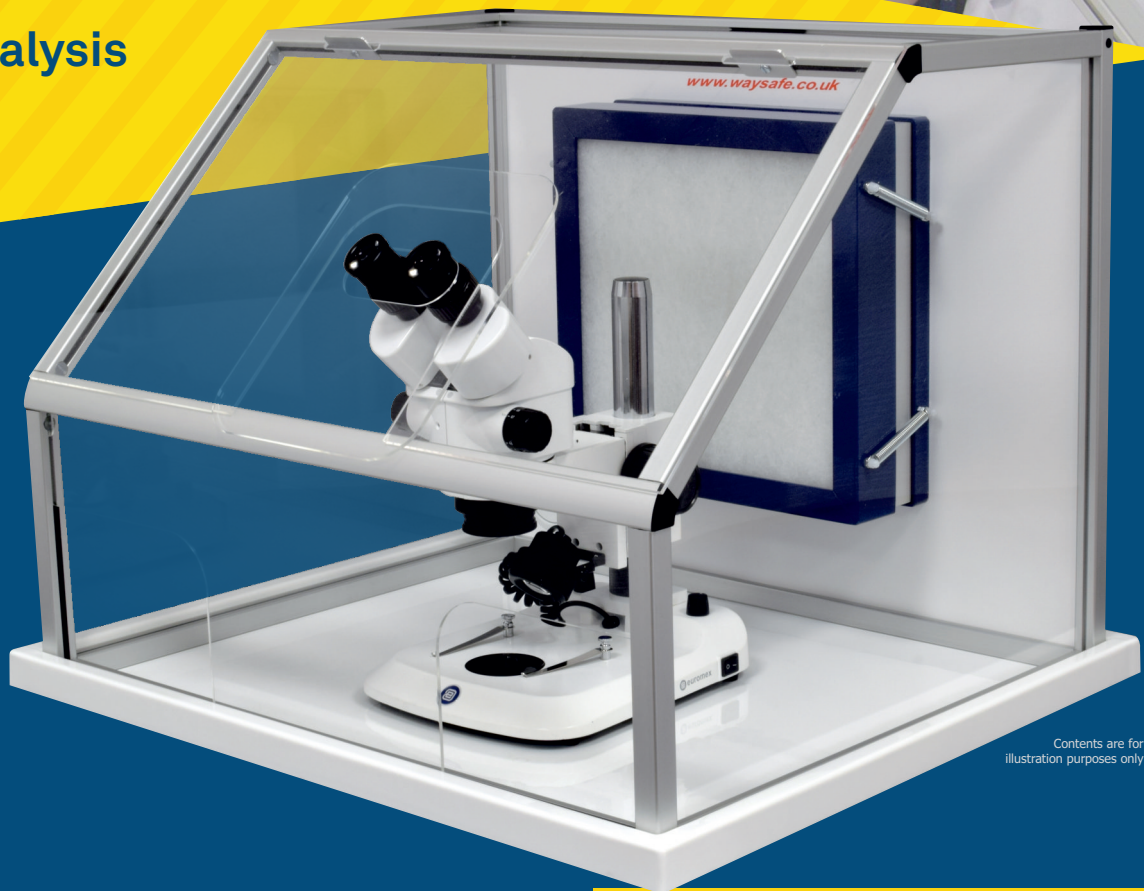

Contents are for illustration purposes only

At 640mm wide the SA640 (previously Waysafe 2) has ample space for most types of microscopes, yet is still small enough to fit on a standard benchtop.

A low noise 2 speed fan draws air through our trusted 300 series H14 HEPA filter before returning it to the lab.



The front panel can be easily removed to provide access the interior for the placing and positioning of equipment. The lower panel has two hand holes for operator access.

Various options available including a range of filters, 120VAC option and an integral composite aluminium base in a range of colours.

- Ductless - air is filtered and clean air circulated back into the workplace
- H14 HEPA filtration (99.999% efficient)
- Easy benchtop installation
- Low maintenance
- 2 speed fan
- Spill tray provided as standard
- Lift off front for ease of access
- Fast global delivery
- Protects operators
- Perfect for microscopy and all types of sample analysis including Asbestos
- Can also be used for other processes media weighing etc.
- Made to order at our factory in the UK

# Technical Specification

## Waysafe SA640

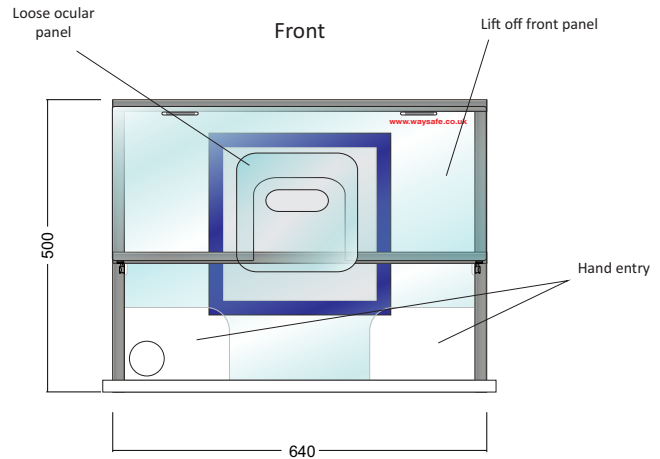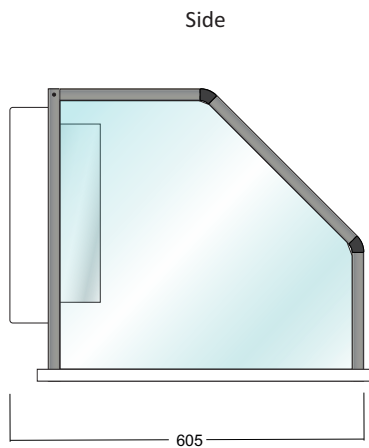

### Technical Specification



<b>Electrical</b>	230VAC 50w, meets LVD and is CE Marked (120VAC available)
<b>Air Speed</b>	High > 0.5m/s, Low 0.3m/s (Face velocity)
<b>External Dim.</b>	Height 500mm, depth 605mm, width 640mm (Tray dimensions 640mm x 680mm)
<b>Filter</b>	Waysafe '300 series' H14 HEPA filter, 99.999% efficiency @ 0.3 micron
<b>Pre filter</b>	Fitted as standard (simple low cost G4 grade filter pad protects the HEPA filter from large particles to extend service life)
<b>Base tray</b>	Push fit moulded spill tray (see options)
<b>Weight</b>	13.9Kg with standard HEPA filter

### Material Specification

<b>Base</b>	High Impact Polystyrene
<b>Glazing</b>	Cast Acrylic Glazing
<b>Frame</b>	Durable Anodised Aluminium
<b>Filter Frame</b>	Powder Coated Aluminium Frame
<b>Ancillaries</b>	Rubber Seals, Durable Hinges

### Optional Extras

<b>Built-in Base</b>	Various colour options including blue, red, white and stainless steel
<b>Left or Right Side Door</b>	Full size opening side door (Only available with Built-in base)
<b>Face Panel</b>	Bespoke panels available to suit varying equipment types
<b>Light</b>	Waysafe TL640 lightbox
<b>Filters</b>	Various filter options include HEPA for particulates and carbons for odours
<b>Variable Speed Fan</b>	Rotary control allows users to steplessly adjust fanspeed for optimum airflow
<b>Hours Run Meter</b>	Electro-mechanical, non-resettable run meter for filter and service monitoring

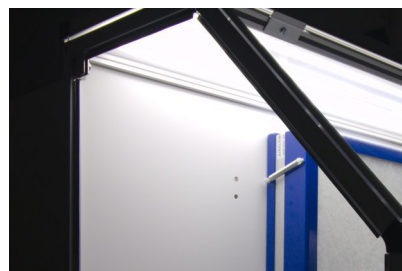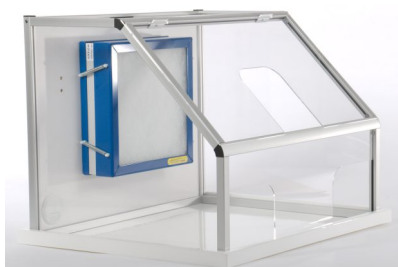

GB2003406





# WolfLabs

**Pricing on any accessories shown can be found by keying the part number into the search box on our website.**

The specifications listed in this brochure are subject to change by the manufacturer and therefore cannot be guaranteed to be correct. If there are aspects of the specification that must be guaranteed, please provide these to our sales team so that details can be confirmed.

**[www.wolflabs.co.uk](http://www.wolflabs.co.uk)**

**Tel : 01759 301142**

**Fax : 01759 301143**

**[sales@wolflabs.co.uk](mailto:sales@wolflabs.co.uk)**

Please contact us if this literature doesn't answer all your questions.